



**Product
Catalogue**







Bathing & Toileting



Aids to Daily Living



Bed & Seating Equipment



Kitchen & Dining

When selecting a Homecraft product, you can expect to experience the difference. Homecraft represents:

Caring spirit

The Homecraft story begins with supporting war veterans, and that heart and compassion extend to all today. Seventy-five years ago, Homecraft identified a need to craft tools for one-of-a-kind situations to keep individuals self-sufficient in their homes. Homecraft is truly the heart of the home.

Dignity

Our intentionally designed products offer Homecraft customers pride and trust to navigate their daily lives. We aim to preserve individuals' confidence and dignity so they may live life fully.

Independence

Homecraft helps individuals declare and craft their own choices and freedom. With an extensive breadth of products, Homecraft has a tool for each situation, giving people the ability to maintain a grasp on life.

Peace of mind

Our attention to detail offers the right solutions to address every need, giving you and your loved ones peace of mind.



Bathing & Toileting



SHOWER CHAIR

Bathe confidently and reduce the risk of falls. The Homecraft Shower Chair allows individuals with limited balance and strength to bathe from a seated, secure position. Included backrest provides upper-body support while bathing. Non-slip feet keep the chair secure during use.

- Maximum user weight of 181.5 kg (400 lb)
- Tool-free assembly for use within minutes
- Height-adjustable from 16 – 25"

Specifications:

Height Adjustable	40.6 x 63.5 cm (16"-25")
Overall height	78.75 x 86.4 cm (31"-34")
Width	40 cm (15.75")
Depth	39 cm (15.35")



**Adjustable legs
with non-slip
feet**



**Drainage holes
prevent water
pooling**



**Handles help
users get up
and down**



Code
7103081

Description
Homecraft Shower Chair



SAVANAH™ SLATTED BATH BOARD

The Savannah Slatted Bath Board has a robust, comfortable, quick-draining design featuring smooth, reinforced plastic slats. Rubber pads under the slats help prevent the bath board from slipping on the bath rim. Brackets can be angled to position the board further into the bath to help secure it.

- **Stainless steel, corrosion-resistant fittings**
- **Brackets include a serrated face and rubber buffers**
- **Slats fully sealed against water ingress**

Code	Length x Width	Weight
091072586	66 x 23.2 cm (26" x 9¼")	2.0 kg (4¼ lb)
091072636	68.6 x 23.2 cm (27" x 9¼")	2.1 kg (4¾ lb)
081124387	71.1 x 23.2 cm (28" x 9¼")	2.2 kg (4½ lb)
081124395	76.2 x 23.2 cm (30" x 9¼")	2.3 kg (5 lb)



190 kg
30 st 419 lb
Maximum user weight

SAVANAH™ SLATTED SHOWER BOARD

The Savannah Slatted Shower Board is similar in function and appearance to the Savannah Slatted Bath Board but supplied with two additional slats for extra comfort and support. It offers a large surface area for users who need to shower on the board and extra width to give users more confidence when transferring. Savannah Slatted Shower Boards may also be used with Savannah Slatted Bath Board Handle.

- **Extra width gives users more confidence when transferring**
- **Stainless steel, corrosion-resistant fittings**
- **Made with high-quality, durable materials**

Code	Length x Width	Weight
091072669	66 x 35.6 cm (26" x 14")	2.8 kg (6 lb)
081124411	68.6 x 35.6 cm (27" x 14")	3.0 kg (6¾ lb)
091072677	71.1 x 35.6 cm (28" x 14")	3.2 kg (7 lb)
081331305	76.2 x 35.6 cm (30" x 14")	3.5 kg (7¾ lb)



190 kg
30 st 419 lb
Maximum user weight

SAVANAH™ SLATTED BATH BOARD HANDLE

The Savannah Slatted Bath Board Handle is an optional accessory for the bath or shower board to ease transfer and offer additional user security. The moulding is covered with a comfortable, soft sleeve to improve handle grip. Fitting the handle may remove the need for a wall-mounted grab rail in many circumstances.



Code	Description	Weight
081124403	Savannah Slatted Bath Board Handle	0.2 kg (7 oz)



LIGHTWEIGHT SUSPENDED BATH SEAT

The Lightweight Suspended Bath Seat is corrosion-resistant and lightweight. With no assembly or adjustments necessary, users can position the bath seat where they need it. Soft caps on the arms protect the edge of the bath and prevent the seat from moving.

- **Lightweight**
- **No assembly required**
- **Constructed from aluminium with a plastic seat**

Specifications:

Seat size	41 x 23 cm (16 x 9")
Overall width	74 cm (29")
Fits baths with internal measurement	58.5-67 cm (23-26 3/8")
Seat height below rim	14 cm (5 1/2")
Weight	1.2 kg (3/4 lb)



Code	Description
091170737	Lightweight Suspended Bath Seat



WHITE LINE SUSPENDED BATH SEAT WITH CUT-OUT

The White Line Suspended Bath Seat with Cut-out features an easy-to-clean plastic seat with a horseshoe cut-out, making personal cleansing easier. Available with or without a backrest.


- **Available as seat only or with backrest**
- **Easy to assemble; supplied flat-packed with tools**
- **Adjustable to suit any bath width**

Specifications:

Seat width	40 cm (16")
Seat depth	40 cm (16")
Seat height below rim	20.5 cm (8")
Overall width adjustment	60-72 cm (23 1/2-28 3/4")
Cut-out	12 x 15 cm (4 3/4 x 6")
Backrest size	47 x 15 cm (18 1/2 x 6")



Code	Description
091477744	Suspended Bath Seat without Back
091356062	Suspended Bath Seat with Back

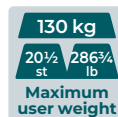


WHITE LINE SWIVELLING BATH SEAT

The White Line Swivelling Bath seat features a lightweight aluminium frame with a swivelling plastic seat and backrest. The seat is designed with a horseshoe front cut-out for personal cleansing and small holes for water drainage. The rotating seat locks at each 90° interval and has side arms to enable easy transfer and provide additional security to the user.

- Aluminium, lightweight frame
- Moulded plastic backrest and seat with cut-out
- Reversible aluminium frame for use on either side of the bath
- Supplied flat-packed with tools

Specifications:

Seat size	41 x 23 cm (16 x 9")
Overall width	74 cm (29")
Fits baths with internal measurement	58.5-67 cm (23-26 3/8")
Seat height below rim	14 cm (5 1/2")
Weight	1.2 kg (3/4 lb)

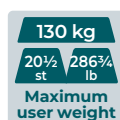

Code	Description
091477769	White Line Swivelling Bath Seat

ALUMINIUM SWIVELLING BATH SEAT

The Aluminium Swivelling Bath Seat is available in three frame widths to suit various bath sizes. The rotating seat locks at each 90° interval and has side arms to enable easy transfer with additional security for the user.

- Lightweight
- Corrosion-resistant
- Seat holes aid water drainage

Specifications:

Seat width	46 cm (18")
Seat depth	39 cm (15 3/4")
Weight	7.2 kg (16 lb)



Code	Frame width
091325836	66.5 cm (26")
091427376	68.5 cm (27")
091325257	71 cm (28")




PADDED TRANSFER BATH BENCH

The Padded Transfer Bath Bench offers an extra-wide bath seat that sits with two legs inside and two legs outside the bath and is designed to support users with limited mobility, helping increase safety during transfers. The user sits securely on the seat and has plenty of room to manoeuvre during a side transfer. The heavily padded seat is comfortable and divided into two sections to allow water to drain.

- Adjustable, aluminium, lightweight frame
- Padded backrest and seat with side handle
- Aluminium frame can be reversed for use on either side of the bath

Specifications:

Adjustable external seat height	55.9-66 cm (22-26")
Maximum clearance from floor	61 cm (24")
Seat size	70 x 40.6 cm (27½ x 16")

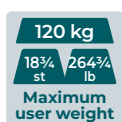

Code	Description
091075662	Padded Transfer Bath Bench



SAVANAH SLATTED BATH SEAT

The Savannah Slatted Bath Seat is a firm, comfortable, quick-draining seat suitable for bathing and showering. Four strong suction cups hold the bath seat securely in place.

- Slats fully sealed against water ingress
- Corrosion-resistant seat
- Available in three heights

Specifications:

Seat size	45.7 x 29 cm (18 x 11¼")
Footprint	34.5 x 30 cm (13½ x 12")



Code	Height	Weight
081124429	15.2 cm (6")	2.2 kg (4¾ lb)
081124437	20.3 cm (8")	2.5 kg (5½ lb)
081124445	30.5 cm (12")	3.2 kg (7 lb)



SAVANAH SHOWER SEAT

The Savannah Shower Seat is a wall-mounted shower seat that conveniently folds neatly against the wall when not used. Made from aluminium and plastic, it is corrosion-resistant and can be easily cleaned. The optional backrest includes arms that fold away independently for side transfers and can be mounted at the desired user height.

- **Generous seat size**
- **Durable**
- **Height-adjustable legs**

Specifications:

Seat width	43 cm (17")
Seat depth	36 cm (14¼")
Seat height	45.5-58.5 cm (18-23")
Distance from wall	43 cm (17")
Distance from wall when folded	13 cm (5")
Weight	3.8 kg (8 lb)



Code	Description
091171511	Savannah Shower Seat
091171529	Optional Back and Arms



PADDED SHOWER SEAT WITH BACK AND ARMS

The Padded Shower Seat with Back and Arms is height-adjustable with a strong, lightweight aluminium frame that makes showering safer for those who cannot stand while bathing. The arms, seat and legs fold quickly away when not used. Includes removable seat and backrest pads.

- **Corrosion-resistant**
- **Lift-up arms**
- **Padded seat, back and arms**

Specifications:

Seat width	46 cm (18")
Seat depth	38 cm (15")
Distance from wall	54.5 cm (21½")
Distance from wall when folded	18 cm (7")
Seat height	46-58 cm (18-22¾")



Code	Frame width
091438589	Padded Shower Seat with Back and Arms



TOOTING PADDED SHOWER SEAT

The Tooting Padded Shower Seat is fitted to the wall and can be folded up when not used. The seat is mounted on an aluminium, stainless steel and plastic frame with adjustable legs fitted with non-slip feet. The legs are adjusted with pin clips and provide extra support and stability. Also available with horseshoe front cut-out.

- Corrosion-resistant frame with padded seat
- Height-adjustable legs with non-slip feet
- Foldable with a latch that holds seat closed

Specifications:

Seat width	39 cm (15½")
Seat depth	33 cm (13")



Code	Description	Weight
081305390	Tooting Padded Shower Seat	4.9 kg (11 lb)
091076884	Tooting Padded Horseshoe Shower Seat	4.8 kg (10½ lb)


TOOTING SLATTED SHOWER SEAT

The White Line Suspended Bath Seat with Cut-out features an easy-to-clean plastic seat with a horseshoe cut-out, making personal cleansing easier. Available with or without a backrest.

- Available as seat only or with backrest
- Easy to assemble; supplied flat-packed with tools
- Adjustable to suit any bath width

Specifications:

Seat width	43.5 cm (17")
Seat depth	29.5 cm (11½")
Seat height adjustment	48.5-59 cm (19-23")
Weight	4.3 kg (9.5 lb)

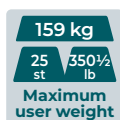

Code	Description
081305382	Tooting Slatted Shower Seat




BUILT-IN HANDLES SHOWER SEAT

The Built-in Handles Shower Seat has fold-down, height-adjustable legs to provide stability. Constructed from easy-to-clean, blow-moulded plastic, the seat has built-in handles to offer the user extra security.

- Legs include non-marking, slip-resistant rubber tips
- Seat includes built-in handles
- Available with steel or aluminium frame

Specifications:

Seat width	49 cm (19¾")
Seat depth	28 cm (11")
Distance from wall	38 cm (15")
Distance from wall when folded	14.5 cm (5¾")

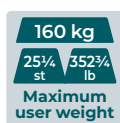

Code	Description	Weight
091171644	Aluminium	4.6 kg (10 lb)
081611029	Steel	5 kg (11 lb)




CONTOUR SHOWER SEAT

The Aluminium Swivelling Bath Seat is available in three frame widths to suit various bath sizes. The rotating seat locks at each 90° interval and has side arms to enable easy transfer with additional security for the user.

- Lightweight
- Corrosion-resistant
- Seat holes aid water drainage

Specifications:

Seat width	46 cm (18")
Seat depth	39 cm (15¼")
Weight	7.2 kg (16 lb)

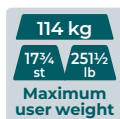

Code	Description
091325265	Contour Shower Seat



WHITE LINE SHOWER SEAT

The height-adjustable White Line Shower Seat has drop-down legs, allowing the seat to be folded away when not used. It is designed for safer independent bathing with a reduced risk of slips and falls. The horseshoe cut-out seat provides easier personal cleansing.

- Aluminium frame with plastic seat
- Horseshoe cut-out seat
- Height-adjustable legs with non-slip feet

Specifications:

Seat width	37 cm (14.5")
Seat depth	40 cm (16")
Seat height adjustment	44-54 cm (17½-21½")
Distance from wall	37 cm (14½")



Code	Description
091477751	White Line Shower Seat

SHERWOOD PLUS BARIATRIC PERCHING AND SHOWER STOOL

The Sherwood Plus Bariatric Perching and Shower Stool is extra-wide and offers generously padded vinyl-covered foam seats for additional comfort. The cross-braced construction and 2.5 cm (1") epoxy-coated steel tubing provide ultra-strong seating support. Choose from two height-adjustable shower stools and two perching stools.

- Height-adjustable legs with pin clips and non-slip feet
- Shower stool available with depth-adjustable backrest
- Perching stools feature a slightly lower frame at the front to give the seat a 5.1-cm (2") incline


Specifications:

Shower stool dimensions:

Seat width	61 cm (24")
Seat depth	43 cm (17")
Height adjustment	50.8-68.6 cm (20-27")
Width between arms	67 cm (26½")

Perching stool dimensions:

Seat width	52.1 cm (20½")
Seat depth	32.4 cm (12¾")
Height adjustment	55.9-73.7 cm (22-29")
Width between arms	61 cm (24")
Footprint at highest setting	77 x 66 cm

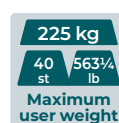

Code	Description
081345206	Shower Stool
081345214	Shower Stool with Adjustable Padded Back

Code	Description
081345180	Perching Stool with Arms and Back
081345198	Perching Stool with Arms and Padded Back



SHERWOOD FOLDING PERCHING STOOL

The Sherwood Folding Perching Stool has a folding frame, making it easy to store when unused. It is height-adjustable and constructed from steel with cushioned vinyl foam seat.

- Eight-adjustable legs with pin clips and non-slip feet
- Available with or without a backrest
- Steel frame with vinyl-covered seat

Specifications:

Seat width	47 cm (18½")
Seat depth	33 cm (13")
Height adjustment	53.5-62.4 cm (21-24½")
Backrest size	30 x 13 cm (11¾")



150 kg
23½ st 330¾ lb
Maximum user weight



Code	Description	Weight
091310846	Standard Stool	8.7 kg (19 lb)
091310853	Standard Stool with Back	10 kg (22 lb)

HEIGHT-ADJUSTABLE PERCHING SHOWER STOOL

The Height-adjustable Perching Shower Stool is available in various styles with or without armrests and backs. The robust frame is manufactured from powder-coated steel, and the grey vinyl seat and back upholstery are fully sealed. The frame is slightly lower at the front to give the seat a slight incline and to enable the user to raise and lower themselves from the stool easily.

- Four different variations of backrests and armrests
- Steel frame with sealed vinyl seats and backrests
- Height-adjustable legs with pin clips and non-slip feet

Specifications:

Seat width	50 cm (20")
Seat depth	28.5 cm (11¼")
Seat height at front	54-70 cm (21¼-27½")
Width between armrests at bottom	46 cm (18")
Width between armrests at top	50 cm (20")

190 kg
30 st 419 lb
Maximum user weight

Code	Description	Weight
091357219	Stool	4.5 kg (10 lb)
091357227	With Arms	5.5 kg (12 lb)
091357235	With Arms and Back	6.9 kg (15 lb)
091357243	With Arms and Padded Back	7 kg (15½ lb)



HEIGHT-ADJUSTABLE BARIATRIC SHOWER STOOL

The Height-adjustable Bariatric Shower Stool has a lightweight aluminium frame and a clip-on plastic seat that can be easily removed for cleaning. The side handles provide additional support when raising or lowering. Optional seat pad available.

- Aluminium frame with vented plastic seat
- Height-adjustable legs with non-slip feet
- Available assembled or flat-packed

Specifications:

Seat width	37 cm (14½")
Seat depth	27 cm (10½")
Seat height	46-58.5 cm (18-23")
Width between arms	50 cm (19¾")
External width	55 cm (21¾")
External depth	37.5 cm (14¾")
Weight	2.8 kg (6 lb)

Code	Description
091325356	Fixed
091564749	Flat-packed
091564756	Optional Seat Pad



Optional seat pad




HEIGHT-ADJUSTABLE ROUND SHOWER STOOL

The Height-adjustable Round Shower Stool features a moulded plastic seat with a sturdy but lightweight frame. The adjustable aluminium legs are finished with rubber ferrules for stability.

- Aluminium frame with round plastic seat
- Height-adjustable legs with non-slip feet; available with two height options
- Available assembled or flat-packed

Specifications:

Seat diameter	34.3 cm (13½")
---------------	----------------

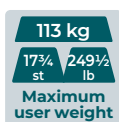

Code	Height	Weight
Assembled		
091075647	40.6-50.8 cm (16-20")	1.5 kg (3¼ lb)
081135342	50.8-61 cm (20-24")	1.7 kg (3¾ lb)
Flat-packed		
091553247	40.6-50.8 cm (16-20")	1.5 kg (3¼ lb)
091553254	50.8-61 cm (20-24")	1.7 kg (3¾ lb)



ALUMINIUM SHOWER STOOL

The Aluminium Shower Stool features a moulded plastic seat with built-in handles for additional user support. Height-adjustable legs are finished with non-marking, slip-resistant rubber tip.

- **Lightweight and corrosion-resistant frame**
- **Vented seat aids water drainage**
- **Available with or without a backrest**

Specifications:

Seat height	38-48 cm (15-19")
Seat width	48.5 cm (19")



091357144



091357151



Code	Description	Weight
091357144	Stool with Backrest	3.1 kg (7 lb)
091357151	Stool	2.5 kg (5½ lb)

WHITE LINE ALUMINIUM CORNER STOOL

The White Line Aluminium Corner Stool features a moulded plastic seat and a lightweight, corrosion-resistant frame. The tripod-shaped frame fits neatly into the corner of a shower cubicle, making optimum use of available space.

- **Lightweight and corrosion-resistant frame**
- **Flat-packed**
- **Assemble with supplied tools**

Specifications:

Seat width	52.5 cm (20½")
Seat depth	38 cm (15")
Seat height	44-54 cm (17¼-21¼")
Footprint	62 x 52 x 52 cm (24½ x 20½ x 20½")
Weight	3.2 kg (7 lb)



Code	Description
091325463	White Line Aluminium Corner Stool



WHITE LINE SHOWER STOOL/CHAIR

The Homecraft White Line Shower Stool/Chair has height-adjustable aluminium frames with moulded plastic seats featuring a cut-out for easier personal cleansing. With optional backrests, armrests and folding, there is a shower stool design for everyone. The modern design frame is lightweight and corrosion-resistant. Supplied flat-packed, the stools are easy to assemble with the tools supplied.

- Lightweight, aluminum frame with moulded plastic seat
- Height-adjustable legs with non-slip feet
- Available with or without armrests and back
- Folding option helps with storage and transportation

Specifications:

Seat width	40 cm (15¾")
Seat depth	40 cm (15¾")



091438670



091438688



091438696



091438704

Code	Description	Ext. dimensions (W x D)
091438670	Shower Stool	47 x 41 cm (18½ x 16")
091438688	Shower Stool with Arms	47 x 41 cm (18½ x 16")
091438696	Shower Chair	51 x 41 cm (20 x 16")
091438704	Folding Shower Chair	44 x 42 cm (17¼ x 16½")



LIGHTWEIGHT CUT-OUT SHOWER CHAIR

The Lightweight Cut-out Shower Chair is height-adjustable and has non-slip ferrules for safety. The frame is anodised aluminium with stainless steel fittings, making it corrosion-resistant. It features an aperture in the seat to make personal cleansing easy and to aid drainage. Armrests provide additional security to the user and may be clipped backward to enable side transfer.

- Anodised aluminium frame with stainless steel fittings
- Heavily padded seat and backrest
- Pan included allowing use as a commode

Specifications:

Seat size	40.5 x 40.5 cm (16 x 16")
Footprint	55.9 x 61 cm (22 x 24")
External width	58.5 cm (23")
Seat height adjustment	44.5-54.5 cm (17½-21½")



Code	Description
081107556	Lightweight Cut-out Shower Chair
081107549	Lightweight Padded Cut-out Shower Chair



BATHMASTER™ DELTIS BATH LIFT

Quickly and easily overcome difficulties getting into and out of the bath with Bathmaster Deltis Bath Lift. This bath lift separates into two parts for easy lifting and carrying. Lithium-ion batteries offer improved performance and battery life. The fully waterproof motor is designed for long-term protection against failure or damage. High-quality materials and construction offer the utmost safety and comfort with no servicing required. It is easy to assemble, and no technical knowledge is needed. Includes bath lift, blue covers, hand controller, one-hour recharger, side flaps and two side-flap protectors.

- Ideal for long-term use
- Compact frame suitable for nearly all baths
- Simple to add and remove from the tub
- Easy to clean and corrosion-resistant

Specifications:

Seat size	50 x 37.5 cm (19¾ x 14¾") excluding side flaps 50 x 69.5 cm (19¾ x 27¾") including side flaps
Backrest size	68 x 35 cm (26¾ x 13¾")
Seat height	6.9-45.5 cm (2¾-18")
Footprint	56 x 29.5 cm (22 x 11½")
Seat weight	6.8 kg (15 lb)
Backrest weight	3.7 kg (8¼ lb)
Maximum user weight	140 kg (308 lb)
Average number of lift cycles	9 (based on 80-kg or 176-lb load)
Guarantee	4 years (all parts)

Bathmaster
deltis



Code
091324144
091324151
091324102
091324110

Description
Bathmaster Deltis Complete UK
Bathmaster Deltis Complete EU
Bathmaster Deltis Heavy Duty Complete UK
Bathmaster Deltis Heavy Duty Complete EU

140 kg
22 st 308¾ lb
Maximum user weight
STANDARD

170 kg
26¾ st 374¾ lb
Maximum user weight
HEAVY DUTY



ATLANTIC COMMODE SHOWER CHAIR

The Atlantic Commode Shower Chair is stylish, easily manoeuvrable and easy to clean. This chair is made from durable plastic-coated steel and is ideal where infection control is a concern. The hygienic clip-on aperture seat, backrest and arms can easily be removed for cleaning or stored away. A toilet roll holder is supplied. It is supplied with a square bowl as standard; however, a disposable pan rack is another option, enabling the use of a disposable pan.


- Upholstery and braked castors are colour-coordinated in blue
- Multifunction armrests are swing back, drop down, detachable and come complete with an easy-to-operate hygienic arm catch
- With optional footrests and disposable pan rack



Specifications:

Overall height		Clearance	44 cm (17½")
Overall width	62 cm (24½")	Width between armrests	51 cm (20")
Overall depth	78 cm (31")	Seat depth	43 cm (17")
Seat height – padded	53.5 cm (20¾")	Castors	13 cm (5")

Code	Description
081533660*	51 cm (20") Bariatric with Footrests
091081280	Standard Square Bowl

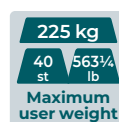

* Available only from US



ALUMINIUM SHOWER COMMUNE CHAIR

The Aluminium Shower Commode Chair is a lightweight, mobile shower chair with four braked castors. The anodised aluminium frame makes it rust-proof and ideal for use in a wet environment. The pushing handle and swivel castors make it easy to manoeuvre, and the blow-moulded plastic seat can be removed for cleaning.

- Anodised aluminium frame with blow-moulded plastic seat
- Padded backrest
- Commode pan included
- Four braked, swivel castors

Specifications:

Seat width	43 cm (17")	Armrest height above seat	17 cm (6¾")
Seat depth	46 cm (18")	External width	56 cm (22")
Seat height	57 cm (22½")	External depth	58.5 cm (23")
Clearance height without pan	49 cm (19¼")	External height	105 cm (41¼")
		Diameter of castors	12.5 cm (5")

Code	Description
091438738	Aluminium Shower Commode Chair



OCEAN WHEELED SHOWER COMMODE CHAIR

The Ocean Wheeled Shower Commode Chair is an attendant-propelled, lightweight mobile shower chair and features a comfortable, contoured horseshoe seat for washing and toileting use and a seat cover with holes for showering. The backrest and armrests are padded for comfort and support. Manufactured from corrosion-resistant materials, it is durable and easy to clean. The aluminium frame can be easily dismantled, ensuring it is compact and lightweight for transport and storage.

- **Lightweight, corrosion-resistant aluminium frame**
- **Two seat options with horseshoe cut-out and full seat**
- **Includes commode pan with lid, footrest, four braked castors and safety belt**



Specifications:

Seat height	52 cm (21")
Seat size	41 x 41 cm (16 x 16")
Size when dismantled	95 x 60 x 24 cm (37½ x 23½ x 9½")
Width between arms	45 cm (18")
Aperture diameter	16 cm (6¼")

Code	Description
081443308*	Ocean Wheeled Shower Commode Chair

* Available only from US



DELUXE SHOWER COMMODE CHAIR

The Deluxe Shower Commode Chair features a lightweight, aluminium frame, making it rust-free and highly durable. Attendant- or self-propelled, this chair offers independence for bathing and toileting. With versatile swing-away, detachable footrests and flip-up arms, it is ideal for use in a shower cubicle or wet room, as a commode or over a standard toilet bowl. The detachable back enables the chair to be partially dismantled for more compact storage.

- Lightweight, aluminium frame with cut-away seat and seat insert
- Heavy-duty version features increased seat width
- Brakes on all four wheels
- Includes square commode pan, footrests and seat insert

Specifications:

Standard chair

Seat height	Fixed height 52 cm (20½")	
	Adjustable height 49.5-57 cm (19½-22½")	
Seat width	45.7 cm (18")	
Seat depth	42.5 cm (16¾")	
Overall height	99 cm (40")	
	Self-propelled	Attendant-propelled
Depth without footrests	84 cm (33")	55 cm (21¾")
Depth with footrests	111 cm (43¾")	81 cm (31¾")
Overall width	69 cm (27")	55.5 cm (21¾")
Aperture size	19 x 36 cm (7½ x 14")	
Front wheel diameter	13 cm (5")	
Rear self-propelling wheels		60 cm (24")
Under seat clearance	Fixed height 46.5 cm (18½")	
	Adjustable height 43-50.5 cm (17-20")	


Heavy-duty chair

Seat height	49.5-57 cm (19½-22½")	
Seat depth	42.5 cm (16¾")	
Overall height	99 cm (40")	
	Self-propelled	Attendant-propelled
Depth without footrests	84 cm (33")	55 cm (21¾")
Depth with footrests	111 cm (43¾")	81 cm (31¾")
Aperture size	19 x 36 cm (7½ x 14")	
Front wheel diameter	13 cm (5")	
Rear self-propelling wheels		60 cm (24")
Under seat clearance	43-50.5 cm (17-20")	

Standard	Description
091183706	Attendant-propelled, fixed height
081135367	Self-propelled, fixed height
091563451	Attendant-propelled, adjustable height

Heavy duty	Description
091563485	Self-propelled
091563477	Attendant-propelled
091563493	Self-propelled


Standard Attendant-propelled

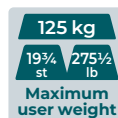

091183706
091563451



Standard Self-propelled

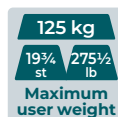

081135367



Heavy-duty Attendant-propelled



091563485



091563477
091563493



Seat width
54 cm (21¼")
54 cm (21¼")
54 cm (21¼")

51 cm (20")
54 cm (21¼")
54 cm (21¼")



DELUXE TILT IN SPACE SHOWER COMMODE CHAIR

The Deluxe Tilt in Space Shower Commode Chair includes all the great features of the Deluxe Shower Commode Chair but has the additional advantages of being height-adjustable and tilting to provide a secure and comfortable position for patients with more complex support requirements. It is a sturdy commode that can tilt the user to provide a secure and comfortable position for washing and toileting. It can be reclined from -5° to 30° with padded neck support that can be positioned at any angle for maximum comfort.

- **Lightweight, aluminium frame with cut-away seat and seat insert**
- **Swing-away, detachable footrests and flip-up arms**
- **Brakes on all four wheels**
- **Includes square commode pan, footrests, seat insert and adjustable neck support**

Specifications:

Seat width	45.7 cm (18")
Seat depth	42.5 cm (16¾")
Seat height	53-60.5 cm (20¾"-23¾")
Overall width	56 cm (22")
Overall height	117 cm (46")
Aperture size	19 x 36 cm (7½ x 14")
Wheel diameter	13 cm (5")
Under seat clearance	43-50.5 cm (17-20")
Weight	10 kg (22 lb)



130 kg
20½ st 286¾ lb
Maximum user weight



Code	Description
091563501	Deluxe Tilt in Space Shower Commode Chair

INFLATABLE SHAMPOO BASIN

The Inflatable Shampoo Basin is designed for comfortable in-bed shampooing with secure, form-fitting support that cushions the head and shoulders. It is lightweight, portable and compact for convenient storage.

- **Built-in inflatable headrest**
- **Made from heavy-duty vinyl**
- **Flexible drainage hose includes stopper for control when emptying**

Specifications:

Width	61 cm (24")
Depth	20 cm (8")
Length	51 cm (20")
Weight	0.68 kg (1½ lb)

Code	Description
091078153	Inflatable Shampoo Basin



INFLATABLE BATH CUSHION

The Inflatable Bath Cushion was designed for comfort and support when bathing. Four suction feet on the underside secure it to the bath to prevent unwanted movement.

Specifications:

Size when inflated	43.2 x 43.2 cm (17 x 17")
Weight	200 g (7 oz)



Code	Description
091077403	Inflatable Bath Cushion

BATH AND SHOWER MAT

The Bath and Shower Mat is a practical solution to bathing security. The basket weave surface reduces the risk of slipping while standing or sitting, and suction cup grips hold the mat in place.

- Made from high-grade rubber
- Multiple drainage holes to drain water quickly
- Available in four sizes

Code	Description	Size
081305408	Regular Bath Mat	57 x 35 cm (22½ x 13¾")
081305416	Large Bath Mat	76 x 35 cm (30 x 13¾")
081107572	Extra Large Bath Mat	57 x 35 cm (22½ x 13¾")
081569755	Regular Shower Mat	54 x 54 cm (21¼ x 21¼")



SURE TREAD BATH AND SHOWER MAT

The non-slip Sure Tread Bath and Shower Mat gives the user total peace of mind. With raised ridges and a unique design, the mat allows the user to stand or sit without slipping on a wet surface. Enjoy the ultimate secure experience of bathing with Sure Tread.

- Made from high-grade medical rubber
- Machine-washable at 40°C (104°F)
- Available in rectangle shape for baths and square for showers

Code	Description	Size
081569821	Bath Mat	78 x 35 cm (30¾ x 13¾")
091077106	Shower Mat	55.5 x 55.5 cm (21¾ x 21¾")



LONG-HANDLED SPONGE

The Long-handled Sponge is mounted on a plastic-coated wire to help reach the legs, feet and back. It is ideal for patients with limited reach.

- Plastic-coated, flexible wire
- Oval foam sponge
- Grey plastic built-up handle



Code	Description
091077585	Length 61 cm (24"), weight 100 g (3.5 oz)
091077569	Length 38 cm (15"), weight 75 g (2.75 oz)

LONG-HANDLED TOE WASHER

The Long-handled Toe Washer has a plastic built-up handle and a flat wire frame covered in a white towelling pad for washing between the toes.

- Plastic-coated, flexible wire
- Grey plastic built-up handle
- Two washing pads included



Code	Description
081305440	Length 71 cm (28"), weight 120 g (4.2 oz)

LONG-HANDLED BATH BRUSH

The Long-handled Bath Brush is evenly balanced and designed to reach the back and legs effectively. The handle shape enables it to be gripped gently without the need for extra pressure. Soft bristles ensure the skin is cleaned without scratching.

- Soft, nylon bristles
- Moulded plastic handle
- Includes hanging loop



Code	Description
091421080	Length 38 cm (15"), weight 150 g (5.3 oz)



FLANNEL BATH STRAP AND WASH MITT

The dual-sided Flannel Bath Strap and Wash Mitt helps the bather reach awkward places. Use the cotton terry surface to wash and the coarse foam reverse side to scour and tone skin. The D-shaped handles create an easy, secure grip, making washing easier.

- Cotton flannel bath strap with D-shaped plastic handles
- Set includes matching hand mitt
- Hand-washable
- Strap dimensions: 80 x 10 cm (31½ x 4"), total length (including handles) 94 cm (37")



Code	Description	Weight
081611078	Flannel Bath Strap and Wash Mitt	100 g (3.5 oz)

DUAL FUNCTION LOTION APPLICATOR

The long-handled Dual Function Lotion Applicator allows users to apply cream and massage their skin simultaneously. Rotating balls ensure even distribution and a pleasant massage. An interchangeable foam head helps when applying thicker creams directly onto the skin. A long, moulded handle makes reaching the back and legs easier. Great for sunscreen application.

- Dual-sided applicator head
- Moulded plastic handle
- No replaceable heads



Code	Description
091077445	Handle length 35 cm (13¾"), weight 150 g (5.3 oz)

SAVANAH MODULAR BATH STEP

The Savannah Modular Bath Step is designed to provide safety and reduce the risk of falling while getting into and out of the bath. The versatile design allows a single step (10-cm or 4-in raise) to be used, or multiple steps can be stacked for additional height. Each step adds 5 cm (2") of height. The large, foam-covered surface gives the user confidence and security.

- Plastic moulded step with textured slip-resistant foam surface
- Modular and stackable (only up to 6 steps should be stacked)
- Available as single step and bulk packed
- Surface area 35.6 x 45.7 cm (14 x 18")
- Footprint 38.1 x 48.3 cm (15 x 19")



Code	Description
081124536	Single Riser, weight 1.3 kg (3 lb)
091162379	Bulk Pack of 12 Steps



BATH GRAB RAIL

The Bath Grab Rail provides a firm handle at the bath side, helping to give confidence and stability when entering and exiting the bath. The coated steel rail is fixed to the side of the bath by a stainless-steel clamp and cushioned on the inside by rubber pads.

- Coated steel frame with stainless steel clamp
- Screw-tight clamp with rubber pads
- Not for use with lipped-style baths

Specifications:

Length	45.7 cm (18")
Height above bath side	34.9 cm (13¾")
Clamp width adjustment	5.5-11.5 cm (2¼-4½")
Weight	2.4 kg (5¼ lb)

Code	Description
081116474	Bath Grab Rail



DELUXE BATH GRAB RAIL

The Deluxe Bath Grab Rail is the perfect solution for safety in the bathroom, offering extra confidence and security while using the bath. The steel-coated rail secures to the side of the tub with a steel clamp and should be checked with each use.

- Coated steel frame with stainless steel clamp
- Screw-tight clamp with rubber pads
- Fits most baths

Specifications:

Length	49.5 cm (19½")
Height above bath side	37.5 cm (14¾")
Clamp width adjustment	7.6-17.8 cm (3-7")
Weight	2.5 kg (5½ lb)

Code	Description
081135433	Deluxe Bath Grab Rail



PADDED BATH GRAB RAIL

The Padded Bath Grab Rail provides added support when entering and exiting the bath or shower. Easy to fit and remove from the tub, the padding offers extra grip ability and comfort for people with limited grip strength. Attractive blue handle is easy to see when using the bath.

- Foam padded, steel frame with adjustable clamp
- Screw-tight clamp with rubber pads
- Fits most baths

Specifications:

Height above bath side	26.7 cm (10½")
Clamp width adjustment	5.7-14 cm (2½-5½")
Weight	2.6 kg (5¾ lb)

Code	Description
081135441*	Padded Bath Grab Rail



* Available only from US



EASYBAR SUCTION GRAB RAIL

The EasyBar Suction Grab Rail offers sturdy and safe support to the user without having to be permanently attached to the wall. The rails assist individuals with sit-to-stand transfers, helping increase safety, balance and stability. The non-slip handle offers greater independence by providing stability when getting in and out of the bath or shower or for general use around the home. Safety lock indicators display green when the rail is securely attached and red when unsecured. The EasyBar was designed for use only on smooth and non-porous surfaces such as tiles, glass, porcelain or fibreglass.

- Equipped with two pressure switches on each end that create a strong vacuum to latch onto surfaces
- Suitable for temporary or permanent use
- Available in two lengths

Code	Description
091325786	35 cm (13¾")
091325802	50 cm (19½")



091325786



091325802


PLASTIC FLUTED GRAB RAIL

The Plastic Fluted Grab Rail is a cost-effective option for bathroom safety. Made from PVC tubes with moulded end fittings, it has circular fixing plates with covering discs that clip into place, hiding the screw heads (screws not included). The tubes are fluted along their length to provide a better grip for wet hands. Available in five sizes.

- White PVC tubes with moulded end fittings and circular fixing plates
- Fluted rail surface
- Available in five lengths

Specifications:

Diameter of tubes	3.8 cm (1½")
Distance from wall	5 cm (2⅓")
Disc diameter	9.1 cm (3⅝")

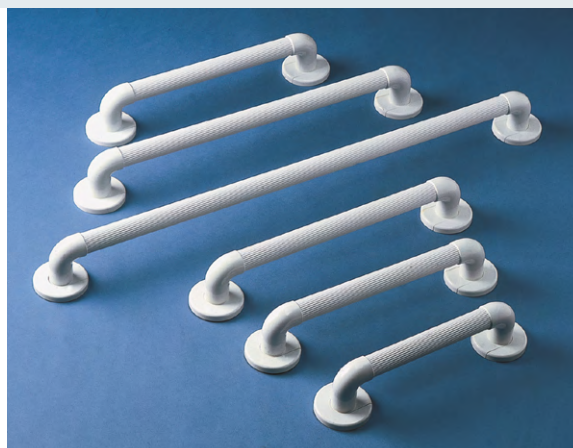

Code	Length	Weight
081306174	30.5 cm (12")	300 g (10.6 oz)
091101641	38.1 cm (15")	370 g (13 oz)
081306190	41 cm (16")	400 g (14 oz)
081306208	45.7 cm (18")	440 g (15.5 oz)
081306216	61 cm (24")	600 g (22.2 oz)



GRIPSURE™ GRAB RAIL

The GripSure Grab Rail is a durable elliptical rail with soft, ribbed mouldings to provide superb grip in wet areas. It is warm to the touch, with a subtle, stylish design suitable for use in various locations around the home. The easy-to-grab ergonomic rail can be mounted horizontally, vertically or diagonally in places to help encourage safety and prevent slips and falls.

- **Ergonomic, elliptical rail with soft, ribbed mouldings**
- **Oval-shaped fixing plates clip into place hiding mounting hardware**
- **Available in four lengths**

Specifications:

Diameter of rails	3.7 cm
Distance from wall	3.7 cm



Code	Length
091187822	30 cm (12")
091187830	40 cm (16")
091188077	45 cm (18")
091188093	60 cm (24")

RINGWOOD GRAB RAIL

The Ringwood Grab Rail is a classic metal rail with a large-diameter tubing for greater comfort and an easier grip. Rail is permanently mounted to the wall and is suitable for indoor and outdoor use.

- **Large diameter: 3.2 cm (1¼") steel tubing**
- **White epoxy or chrome-plated finish**
- **Available in four lengths**



Code	Description	Length
091100874	White Epoxy-coated	30.5 cm (12")
091101138	White Epoxy-coated	45.7 cm (18")
091101195	White Epoxy-coated	61 cm (24")
091101211	White Epoxy-coated	91.4 cm (36")
091101120	Chrome	30.5 cm (12")
091101187	Chrome	45.7 cm (18")
091101203	Chrome	61 cm (24")
091101229	Chrome	91.4 cm (36")



STEEL GRAB BAR

The Steel Powder-coated Grab Bar is made of epoxy steel tubing for durability and strength. It is suitable for indoor and outdoor use. Add a rail near the threshold to help manage stairs or in the bathroom to help prevent slips and falls.

- 2.5 cm (1") steel tubing rail
- White epoxy finish
- Available in two lengths

Specifications:

Width	61 cm (24")
Depth	20 cm (8")
Length	51 cm (20")
Weight	0.68 kg (1½ lb)



Code	Length
091393230	44.5 cm (17½")
091393248	61 cm (24")

CHROME-PLATED STEEL TEXTURED GRAB RAIL

The Chrome-plated Steel Textured Grab Rail is durable and practical with a patterned, indented grip to help prevent hands from slipping. The tubing is a comfortable diameter and stands off from the wall by 3.8 cm (1½").

- 2.5 cm (1") steel tubing rail with textured grip
- Chrome-plated finish
- Available in five lengths



Code	Length
091101427	30.5 cm (12")
081570464	40.6 cm (16")
081570472	45 cm (18")
081570480	61 cm (24")
091101484	81.2 cm (32")



POLISHED STRAIGHT/CURVED STAINLESS STEEL GRAB BAR

The stylish Polished Straight Stainless Steel Grab Bar has been designed to offer stability and support in the bathroom or around the home. The rails are durable and rust-resistant with stainless steel construction, and a polished mirror finish provides a high-quality, modern look to complement any bathroom. Oval-shaped fixing plates include attractive covers to hide unsightly screw heads.

- Rust-resistant, stainless-steel rail with oval fixing plates and screw covers
- Available in straight or curved rail design
- Both designs available in three lengths

Specifications:

Diameter of rails	3.2 cm (1¼")
Distance from wall	5 cm (2")

Code	Description	Length
Straight		
091329101	Straight	30.5 cm (12")
091329119	Straight	45 cm (18")
091329127	Straight	61 cm (24")

Curved		
091329143	Curved	30.5 cm (12")
091329150	Curved	45 cm (18")
091329168	Curved	61 cm (24")



SAVANAH RAISED TOILET SEAT

The Savannah Raised Toilet Seat features a gently contoured surface for extra comfort and slopes down slightly from the back to the front for easier lowering and raising. A large rear cut-out gives clearance for the coccyx and combined with the gentle downward slope of the seat, fosters rear personal hygiene and an expansive front cut-out. The deep, built-in splash guard minimises soiling. Fully sealed, the seat is resistant to stains and odours and is easy to clean.

Two L-shaped moulded plastic brackets have bonded non-slip rubber buffers to keep the toilet firmly in place. The brackets are easily adjusted to fit most toilet shapes and sizes using two large, fluted turning knobs, making fitting and removing the seat quick and easy. No front bracket or front-turning knob to avoid snagging skin or clothes.

- **Lightweight, plastic construction for strength and durability**
- **Fully sealed, resistant to stains and odours**
- **Available in three seat heights**
- **Available with or without lid**

Specifications:

Overall length	40.7 cm (16")
Overall width with brackets	37 cm (14¾")
Aperture	26.6 x 20.7 cm (10½ x 8¼")

Code	Description	Height	Weight
091079185	Without Lid	5 cm (2")	0.8 kg (1¾ lb)
091079276	Without Lid	10 cm (4")	1 kg (2¼ lb)
081305549	Without Lid	15 cm (6")	1.15 kg (2½ lb)
081305515	With Lid	5 cm (2")	1.2 kg (2¾ lb)
081305531	With Lid	10 cm (4")	1.4 kg (3 lb)
081305556	With Lid	15 cm (6")	2 kg (4½ lb)



PADDED RAISED TOILET SEAT

The Padded Raised Toilet Seat features a soft cushion to provide extra comfort and raise the effective seat height. The white vinyl covering protects the cushion is easy to clean and is attachable to the existing toilet seat with four hook-and-loop straps.

- **White, vinyl-covered cushion with hook and loop fasteners**
- **Attaches to existing toilet seat**
- **Available in two seat heights**
- **Available with or without lid**

Specifications:

Width	36 cm (14")	
Length	41 cm (16")	
Code	Height	Weight
091079540	5 cm (2")	270 g (9.5 oz))
091079557	10 cm (4")	380 g (13.4 oz)



RAISED TOILET SEAT WITH ARMRESTS

The Raised Toilet Seat with Armrests fosters increased independence, notably by compensating for impaired movements. The seat with detachable armrests comes with a lid, that can be removed for easy cleaning. The versatile design includes a choice of three fixed seat heights to suit the user's needs.

- Moulded plastic base includes lift-up ABS armrests and a steel frame
- Cross bars adjust in width for each seat position
- Existing toilet seat must be removed prior to installation

Specifications:

Seat width	36 cm (14")
Seat depth	42.5 cm (16 ² / ₃ ")
Overall width	50 cm (19 ¹ / ₂ ")
Overall depth	47.5 cm 18 ² / ₃ ")
Overall height	38 cm (15")
Seat height adjustments	6 cm (2"), 10 cm (4") and 15 cm (6")
Aperture size	21 x 27.5 cm (8 ¹ / ₄ x 10 ² / ₃ ")

Code	Description
081703990	Raised Toilet Seat with Armrests



136 kg
21¹/₄ st 300 lb
Maximum user weight

RAISED TOILET SEAT WITH BUILT-IN ARMS

The Raised Toilet Seat with Built-in Arms features a moulded plastic seat with removable aluminium arms and a large front cut-out for easy personal cleansing. Easily fitted to the toilet bowl with three adjustable brackets. The arms are covered with soft foam padding for comfort and grip.

- Moulded plastic, raised toilet seat with built-in arms
- Removable, foam-padded arms
- Three adjustable brackets attach seat to toilet

Specifications:

Height	10 cm (4")
Aperture size	21 x 25.5 cm (8 ³ / ₄ x 10")
Width between armrests	44.5 cm (17 ¹ / ₂ ")
Height of armrest above seat	19 cm (7 ¹ / ₂ ")
Weight	7.5 kg (16 ¹ / ₂ lb)



Code	Description
091079722	Raised Toilet Seat with Built-in Arms

180 kg
28³/₄ st 397 lb
Maximum user weight



FOLDAWAY GRAB RAIL

The Foldaway Grab Rail is made from powder-coated steel and designed for use next to a toilet to increase safety while lowering and raising. The wall bracket attaches permanently to the wall and supports the rail in an upright position when unused. It is available with a height-adjustable supporting leg, which provides additional user weight capacity and features a non-marking rubber tip and locking mechanism when the rail is in an upright position.

- **White, powder-coated steel rail**
- **Attaches to wall with steel wall bracket**
- **Available with and without supporting leg**

Specifications:

Fold Away Grab Rail

Dimensions of the wall bracket 28.2 x 10.4 cm (11 x 4")

Rail length 75 cm (29½")

Fold Away Grab Rail with Supporting Leg

Height adjustment 80.4-112.8 cm (31½-44½")

Dimensions of the wall bracket 20.7 x 11 cm (8¼ x 4½")

114 kg
17¾ st 251½ lb
Maximum
user weight

125 kg
19¾ st 275½ lb
Maximum
user weight

Code	Description
091357029	Fold Away Grab Rail
091357037	Fold Away Grab Rail with Supporting Leg



DEVON FOLDING SUPPORT RAIL

The Devon Folding Support Rail is a sturdy folding rail that projects at right angles to the wall. The rail is held in the upright position by a spring-loaded cam and is easily folded down.

- **White, powder-coated steel rail**
- **Attaches to wall with steel wall bracket**
- **Available in standard and short length**

Specifications:

Rail diameter 2.5 cm (1")

Backplate dimension 30 x 12.5 cm (11¾ x 5")



67 kg
10½ st 147¾ lb
Maximum
user weight

Code	Description	
081305465	Standard — 76 cm (30")	3.0 kg (6½ lb)
081305457	Short — 55 cm (21½")	2.6 kg (5¾ lb)



DEVON DELUXE FOLDING SUPPORT RAIL

The Devon Deluxe Folding Support Rail offers the same features as the Devon Folding Support Rail but with a larger diameter rail for easier gripping, a spring mechanism to hold it upright, as well as a locking grip. Available as a standard folding rail with a supporting leg or as a floor-mounted rail, this product is available with multiple options to meet the user's specific requirements.

84 kg
13½ st 185½ lb
Maximum user weight

- **White, powder-coated steel rail with larger diameter rail**
- **Spring mechanism and lock hold rail in the upright position**

Specifications:

Rail diameter 3.2 cm (1¼")
 Backplate dimension 30 x 12.7 cm (11¾ x 5")

Code	Rail length	Weight
081305473	Standard — 77 cm (30½")	4.1 kg (9 lb)
091078567	Short — 56 cm (22")	3.4 kg (7½ lb)



DEVON DELUXE FLOOR-MOUNTED FOLDING RAIL

The Devon Deluxe Floor-mounted Folding Rail has a large diameter rail for easier gripping and a spring mechanism to hold it in the up position. It is mounted on a floor-fixing post for positioning versatility with a clamp to adjust by height.

127 kg
20 st 280 lb
Maximum user weight

Specifications:

Floor post height 111 cm (44")
 Post height adjustment 58.4-101.2 cm (23-40")

Code	Rail length	Weight
081305499	Standard — 76.5 cm (30")	8.4 kg (18½ lb)
091078617	Short — 55.9 cm (22")	6.8 kg (15 lb)



DEVON DELUXE FOLDING SUPPORT RAIL WITH LEG

The Devon Deluxe Folding Support Rail with Leg features a support leg to offer strength and stability during lowering and raising, lessening the strain on the wall. The height-adjustable leg includes a non-slip ferrule.

Specifications:

Height adjustment 90-112.9 cm (35½-44½")
 Backplate 30 x 12.7 cm (11¾ x 5")

Code	Rail length	Weight
081305481	77.5 cm (30½")	5 kg (11 lb)



67 kg
10½ st 147¾ lb
Maximum user weight



DEVON ELITE FOLDING SUPPORT RAIL

The Devon Elite Folding Support Rails were designed to offer superb support, providing the user with greater confidence, security and independence in the bathroom. Manufactured from light yet strong aluminium, the elliptical arms give the user the ultimate comfort and grip while lowering to or raising from the toilet.

- Lightweight, aluminium frame with elliptical-shaped rail
- Arm-locking mechanism and spring mechanism support weight of the rails as lowered
- Available in three versions: standard folding rail, rail with supporting leg or floor-mounted rail

Specifications:

Wall-mounted Folding Support Rail

Rail width	5 cm (2")
Rail depth	2.8 cm (1")
Back plate	23 x 12 cm (9 x 4½")
Projection when folded	17.5 cm (7")

127 kg
20 st 280 lb
Maximum user weight

Floor-mounted Folding Support Rail

Rail width	5 cm (2")
Rail depth	2.8 cm (1")
Height from floor to top of arm	62.5-98.5 cm (24½-38¾")

127 kg
20 st 280 lb
Maximum user weight



081446368



081437144



091173145

Code	Description	Rail length
081446368	Wall-mount Folding Support Rail	64.5 cm (25½")
081437144	Floor-mounted Folding Support Rail	64.5 cm (25½")
091173145	Folding Support Rail with Adjustable Leg	64.5 cm (25½")

RINGWOOD WALL-TO-FLOOR GRAB RAIL

The Ringwood Wall-to-Floor Grab Rail is standard floor-mounted grab rail with larger-diameter steel tubing for greater control and easier grip. A white epoxy coating makes it suitable for indoor or outdoor use. Rail fixes to the floor and wall measuring 76 cm (30") from the floor and 78 cm (30¾") from the wall.

- 3.2 cm (1¼") white, epoxy-coated steel with large rail diameter
- Wall- and floor-mounted
- Suitable for indoor and outdoor use

Code	Description
081330901	Ringwood Wall-to-Floor Grab Rail



STIRLING TOILET FRAME

The Stirling Toilet Frame features a sturdy, powder-coated metal frame with a clip-on, one-piece moulded toilet seat that is slightly angled forward to make lowering and raising easier. The contoured seat has a partial front cut-away to facilitate personal hygiene. The width-adjustable option features a pin-clip system, enabling the frame to fit a variety of toilet styles and be completely assembled or disassembled without tools. The deluxe version includes padded armrests on the side rails for additional comfort and easier grip.

- Frame legs feature non-slip rubber feet
- Pin-clip system for height and width adjustments leaves no protrusions to snag skin or clothing
- Flat-packed, assembly required
- Available in three versions: standard, width-adjustable or width-adjustable with padded arms



Code	Description
091080027	Standard Height-adjustable Frame
091080068	Standard Height-/ Width-adjustable Frame
091080076	Deluxe Height-/ Width-adjustable Frame with Padded Arms



TOILET SAFETY FRAME

Frame that provides a secure hand grip and stability for users. The frame is secured to the toilet using the existing seat, so no additional parts or locking mechanisms are necessary. The height-adjustable legs suit various toilet designs.

- Corrosion-resistant aluminium frame
- Height-adjustable legs with non-slip feet
- Moulded armrests

Specifications:

Height of arms	64-75 cm (25-29½")
Width between legs	40-55 cm (16-22")
Width between arms centre to centre	53-63 cm (21-25½")
Overall depth	47 cm (18½")

Code	Description
091078872	Toilet Safety Frame



UNI-FRAME FOLDING TOILET SURROUND

The Uni-Frame Folding Toilet Surround, like the Uni-Frame Folding Toilet Frame, has a folding frame design but without the toilet seat. The support rail is height-adjustable and can be folded flat for storage or transport. It offers convenient support for additional bathroom safety and security.

- Coated metal frame, folds flat for storage
- Height-adjustable legs with non-slip feet
- Padded backrest

Specifications:

Height	62-74.5 cm (24½-29¾")
Width	50 cm (19½")
Depth	42 cm (16½")
Weight	5.9 kg (13 lb)

Code	Description
081124551	Uni-Frame Folding Toilet Surround



HEIGHT-ADJUSTABLE TOILET SURROUND

The Height-adjustable Toilet Surround, available in steel or aluminium, provides extra support and safety when using the toilet. The height-adjustable frame features moulded plastic armrests for additional comfort.

- Available with steel or lightweight aluminium frame
- Height-adjustable legs with non-slip feet
- Moulded plastic armrests

Specifications:

Steel

Height	69.5-85 cm (27½-22½")
Width	59 cm (23¼")
Depth	46 cm (18")
Weight	3 kg (6¾ lb)

160 kg
25½ st 352¾ lb
Maximum user weight

Aluminium

Height	66-76 cm (26-30")
Width	60.5 cm (23¾")
Depth	47 cm (18½")
Weight	2.1 kg (4½ lb)

125 kg
19¾ st 275½ lb
Maximum user weight

Code

091325711

Description

Height-adjustable Toilet Surround — Steel

091078880

Height-adjustable Toilet Surround — Aluminium



FULLY ADJUSTABLE TOILET SURROUND

The Fully Adjustable Toilet Surround, available with a steel or aluminium frame, is both width- and height-adjustable, making it suitable for a variety of toilet styles and users' needs. Padded handgrips provide additional comfort and support.

- Available with steel or lightweight aluminium width-adjustable frames
- Height-adjustable legs with non-slip feet
- Padded armrests
- Supplied flat-packed

Specifications:

Steel

Height	62.5-78.8 cm (25-31")
Width	51-56 cm (20-22")
Depth	49.5 cm (19½")
Weight	5 kg (11 lb)

190 kg
30 st 419 lb
Maximum user weight

Aluminium

Height	62.5-78.8 cm (25-31")
Width	51-56 cm (20-22")
Depth	49.5 cm (19½")
Weight	2.8 kg (6¼ lb)

190 kg
30 st 419 lb
Maximum user weight

Code

081124577

Description

Fully Adjustable Toilet Surround — Steel

081534684

Fully Adjustable Toilet Surround — Aluminium





091080407



091325604



091309798



091438993

STACKING COMMODE

The Stacking Commode is a sturdy, metal frame commode with a hinged black plastic toilet seat over a commode pan. Available as a fixed-height or height-adjustable commode, it features a vinyl-covered padded lift-off seat and padded backrest. It can be conveniently stacked if needed.

- **Black, metal frame**
- **Hinged black plastic toilet seat**
- **Height-adjustable version features pin clips and plastic feet**
- **Includes commode pan and lift-off seat**

Specifications:

Stacking Commode

Width at base	58 cm (23 ³ / ₄ "
Depth at base	50 cm (19 ³ / ₄ "
Seat height	47 cm (18 ¹ / ₂ "
Aperture size	24 x 24 cm (9 ¹ / ₂ x 9 ¹ / ₂ "
Weight	6.7 kg (14 ³ / ₄ lb)



Adjustable Stacking Commode

Seat depth	43 cm (17"
Seat height (aperture seat)	48-58 cm (19-23"
Seat height (padded)	51.5-61.5 cm (20 ¹ / ₄ -24 ¹ / ₄ "
Width between arms	47 cm (18"
Overall width	53 cm (21"
Overall depth	45.5 cm (18"
Overall height	80-90 cm (31 ¹ / ₂ -35 ¹ / ₂ "
Aperture size	24 x 24 cm (9 ¹ / ₂ x 9 ¹ / ₂ "



HEIGHT-ADJUSTABLE COMFORT COMMODE

The Height-adjustable Comfort Commode with a sturdy metal frame, non-slip feet, padded backrest and padded seat cover offers comfort and stability within any home or institution. The padded blue seat cover and clip-on black toilet seat are easily removed for cleaning. The deluxe version has a padded blue toilet seat that provides the option of additional comfort. The commode pan with a handle is easy to remove for waste disposal. Also available with removable armrests.

- **Blue, metal frame**
- **Clip-on black plastic toilet seat**
- **Height-adjustable legs feature pin clips and non-slip feet**
- **Deluxe version includes padded armrests and padded toilet seat**
- **Includes commode pant**

Specifications:

dth	45 cm (17 ³ / ₄ "
Seat depth	43 cm (17"
Seat height (padded)	48.5-63.5 cm (19-25"
Seat height (classic)	45-60 cm (17 ³ / ₄ -23 ¹ / ₂ "
Height to top of armrest	59-79 cm (23 ¹ / ₄ -31"
Overall width	55 cm (21 ¹ / ₂ "
Overall depth	55 cm (21 ¹ / ₂ "
Overall height	83-98 cm (32 ³ / ₄ -38 ¹ / ₂ "
Commode pan	33 x 16 cm (13 x 6 ¹ / ₂ "



Standard	Description
091080407	Stacking Commode
091325604	Adjustable Stacking Commode

Standard	Description
091309798	Standard
091325828	Standard with Removable Armrests
091438993	Deluxe with Padded Toilet Seat and Armrests
091325760	Spare Commode Pan



DELUXE COMFORT FOLDING COMMODE

The Deluxe Comfort Folding Commode is made from epoxy-coated steel and conveniently folds for storage and transportation. The backrest and toilet seat are padded to provide additional comfort for the user.

- Epoxy-coated steel frame with plastic feet
- Padded toilet seat, backrest and seat cover
- Includes removable, padded seat cover, commode pan and commode pan lid

Specifications:

Seat width	46 cm (18")
Seat depth	45 cm (17¾")
Seat height	48 cm (19")
External width	51 cm (20")
External height	90 cm (35½")
Width between arms	46 cm (18")
Weight	7 kg (15½ lb)

Code	Description
091081017	Deluxe Comfort Folding Commode



FOLDING COMMODE AND TOILET SURROUND

The Folding Commode and Toilet Surround is a lightweight aluminium surround with a clip-on toilet seat that can be easily removed for cleaning. The toilet surround with padded arms provides extra support and safety when using the toilet and eliminates the need for a grab rail. It conveniently folds for storage or transport.

- Lightweight, aluminium frame
- Plastic, clip-on toilet seat
- Height-adjustable legs with non-slip feet
- Includes removable toilet seat and commode pan

Specifications:

Seat width	34.5 cm (14")
Seat depth	37 cm (14½")
Seat height adjustment	46-56 cm (18-22")
Width between arms	46 cm (18")
Overall width	55 cm (21¾")
Overall depth	49 cm (19¾")

Code	Description
091439348	Folding Commode and Toilet Surround




CHROME-PLATED WHEELED COMMODE

The Chrome-plated wheeled commode features a standard steel frame with detachable armrests for easy side transfers and a detachable backrest for storage. Push handles and swivelling castors allow the chair to be easily manoeuvred over the toilet. The blow moulded aperture seat is hygienic and easy to clean.

- Blow-moulded aperture seat with removable padded cover
- Detachable armrests with moulded supports
- Available with or without footrests

Specifications:

Seat width	43 cm (17")
Seat depth	45 cm (17¾")
Seat height (padded)	54 cm (21¼")
Seat height (aperture)	50 cm (19¾")
Overall height	87.5 cm (34½")
Overall width	52 cm (20½")
Overall depth	57 cm (22½")
Top of padded seat to armrest	22 cm (8½")
Floor to top of armrest	75 cm (29½")
Diameter of castors	12 cm (4¾")
Weight	13.5 kg (29¾ lb)

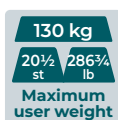

091357094



091325729



Code	Description
091357094	Commode Round Bowl without Footrests
091325729	Commode with Footrests
091439827	Optional Footrests (Pair)



HEIGHT-ADJUSTABLE WHEELED COMMODE



The Height-adjustable Wheeled Commode is manufactured from durable chrome-plated steel and is available in three seat widths. It is height-adjustable and includes four swivel brake castors and a push handle to allow the commode to be wheeled over a toilet if needed. Side transfers are easier and safer with swinging drop arms. The easy-to-clean, hygienic, blow-moulded aperture seat includes a removable padded seat cover so the commode can be used as a transfer chair.

- Chrome-plated steel frame with swivelling, braked castors, and push handle
- Detachable armrests with moulded supports, and detachable back
- Three seat width options
- Includes square commode pan and hook-on footrests

Specifications:	18" mobile	20" mobile	22" mobile
Seat width	43.5 cm (17")	43.5 cm (17")	43.5 cm (17")
Seat depth	46 cm (18")	46 cm (18")	46 cm (18")
Seat height (padded)	54.5-64.5 cm (21½-25½")	54.5-64.5 cm (21½-25½")	54.5-64.5 cm (21½-25½")
Seat height (aperture)	50-60 cm (19½-23½")	50-60 cm (19½-23½")	50-60 cm (19½-23½")
White aperture seat orifice size	28.5 x 18 cm (11¼ x 7")	28.5 x 18 cm (11¼ x 7")	28.5 x 18 cm (11¼ x 7")
Width between arms	46 cm (18")	50 cm (20")	56 cm (22")
Top of padded seat to armrest	22 cm (8.5")	22 cm (8.5")	22 cm (8.5")
Floor to top of armrest	76-86 cm (30-34")	76-86 cm (30-34")	76-86 cm (30-34")
Overall width	54.5 cm (21½")	60 cm (23½")	66 cm (26")
Overall depth	62 cm (24½")	62 cm (24½")	62 cm (24½")
Overall height	95-105 cm (37½-41½")	95-105 cm (37½-41½")	95-105 cm (37½-41½")
Wheels	12 cm (5")	12 cm (5")	12 cm (5")
Weight	17.5 kg (38½ lb)	18.5 kg (40¾ lb)	19 kg (42 lb)

Standard	Description
091081124	18" Seat width
091438936	20" Seat width
091438134	22" Seat width



FOLDING BOTTOM WIPER

The Folding Bottom Wiper features a moulded plastic head with recessed serrations to grip paper. It folds into three sections for easier storage. Designed to be used from the front.

Specifications:

Length	26.5 cm (10½")
Fold-down length	12 cm (4¾")
Weight	0.8 kg (1¾ lb)

Code	Description
081305622	Folding Bottom Wiper



BOTTOM WIPER

Designed by the Bath Institute for Medical Engineering, the Bottom Wiper has a moulded plastic head with recessed serrations to grip paper. It is mounted on a long, lightweight, curved handle that assists in reaching areas where wiping is difficult.

Specifications:

Length	26.5 cm (10½")
Weight	0.7 kg (1½ lb)

Code	Description
081305614	Bottom Wiper



SAVANAH SPLASH GUARD

Simple to fit, the Savannah Splash Guard clips onto a toilet bowl, a raised toilet seat or a toilet frame to prevent embarrassing and unhygienic accidents. Simple to remove for cleaning and getting on and off the toilet.

Code	Description
081441377	Savanah Splash Guard



GENERAL PURPOSE BEDPAN

The General Purpose Bedpan is made from polypropylene and can be cleaned with steam, hot water or chemicals and has a 1-litre capacity.

- **Lightweight, aluminium frame with cut-away seat and seat insert**
- **Swing-away, detachable footrests and flip-up arms**
- **Brakes on all four wheels**
- **Includes square commode pan, footrests, seat insert and adjustable neck support**

Specifications:

Width	30 cm (11¾")
Length	38 cm (15")
Depth	10 cm (4")
Weight	1.1 kg (2½ lb)



Code	Description
091081876	General Purpose Bedpan

ECONOMY MALE OR FEMALE URINAL

The Economy Male Urinal is moulded from transparent plastic and features a snap-on lid and handle. The long neck ensures it is easy to use. It has a 1-litre capacity.

The Economy Female Urinal has a handle for ease of use. It is anatomically designed to use when sitting, standing or lying down. It has an 800-ml capacity with graduations for output.



081610872



091082049

Code	Description	Weight
081610872	Male Urinal	0.1 kg (3.5 oz)
091082049	Female Urinal	0.1 kg (3.5 oz)



Aids to daily living



DRESSING STICK

The Dressing Stick helps individuals put on and take off garments they cannot easily reach or that require movements that are difficult, such as socks, jumpers and jackets. Each dressing stick is equipped with a rubber tip on one end and a double-wire hook on the other.

Specifications:

Length	52 cm (21")
Weight	36 g (1.3 oz)

Code	Description
091095686	Dressing Stick



KINGS BUTTON HOOK

The Kings Button Hook is a shaped stainless-steel wire constructed to assist with one-handed dressing. Simply pass the hook through a buttonhole and over the button, then pull the button back through the buttonhole.

Specifications:

Wire length	9.7 cm (3¾")
Handle	12 cm (4¾")
Weight	75 g (2.6 oz)



Code	Description
091095678	Kings Button Hook



BUTTON HOOK AND ZIP PULLER

The Button Hook and Zip Puller is a 2-in-1 tool that assists individuals in buttoning clothing and pulling zippers. A button hook on one end and a zip puller on the other help with various dressing tasks.

Specifications:

Length	24.5 cm (9¾")
Handle width	1.5 cm (¾")
Wood length	17 cm (6¾")
Weight	8 g (0.3 oz)



Code	Description
091095645	Button Hook and Zip Puller



PLASTIC SHOEHORN

The Plastic Shoehorn features a long, curved handle designed to extend reach when putting on shoes. Includes an attached loop for storage.

Specifications:

Length	43 cm (17")
Weight	37 g (1.3 oz)



Code	Description
091095140	Plastic Shoehorn



EXTRA LONG PLASTIC SHOEHORN

The Extra Long Plastic Shoehorn offers additional reach, eliminating bending and stooping when putting on footwear. Includes an attached loop for storage.

Specifications:

Length	56 cm (22")
Maximum width	4.5 cm (1¾")
Weight	45 g (1.6 oz)



Code	Description
091318070	Extra Long Plastic Shoehorn



CHROME-PLATED METAL SHOEHORN

The Chrome-plated Metal Shoehorn features a chrome metal finish on a long-stemmed shoehorn. The blade is mounted on a spring to allow it to bend in relation to the handle when putting on shoes. The plastic-moulded handle has a comfortable grip, making it easy for the user to hold.

Specifications:

Length	60 cm (23¾")
Weight	115 g (4 oz)

Code	Description
091095199	Chrome-plated Metal Shoehorn




COATED METAL SHOEHORN

The Coated Metal Shoehorn features a sturdy, long design with epoxy-coated steel and a soft, hook-shaped PVC handle.

Specifications:

Length	60 cm (23¾")
Weight	233 g (8.2 oz)

Code	Description
091095173	Coated Metal Shoehorn

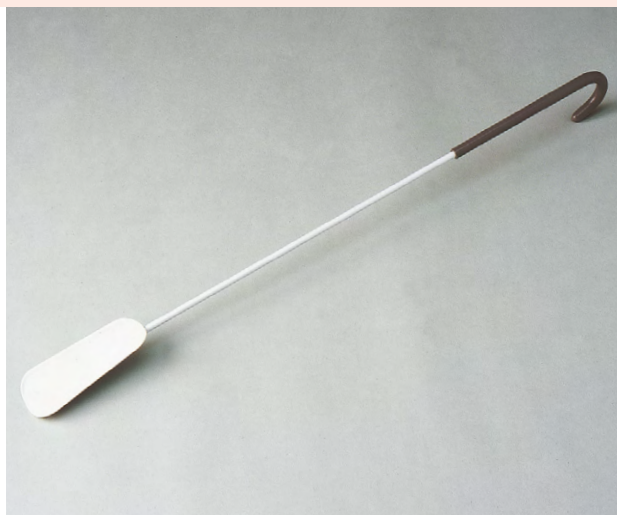

DORKING STOCKING AID — STANDARD OR SIDE SLOT


The Dorking Stocking Aid is a simple, inexpensive device that assists the user when putting on stockings or socks, especially when bending to reach the foot is difficult. It features a flexible plastic gutter with long cotton tapes attached to the top. Feed the stocking or sock through the plastic and hold it in place with the side notches. Place the foot in the open end of the stocking and pull it over the foot and up the leg using the tapes. Available with or without side slots.

Specifications:

Cord length	92 cm (36")
Gutter length	22.9 cm (9")
Weight	55 g (2 oz)



Code	Description
091095421	Standard (Non-Retail)
081305705	Standard (Retail)
091095454	Side Slot (Non-Retail)

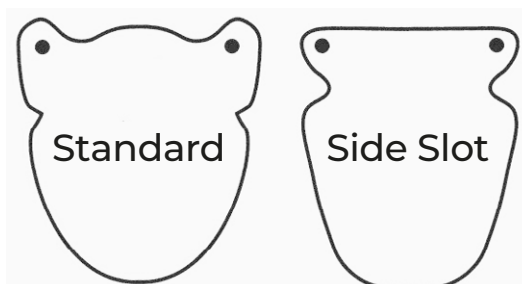

BREVETTI SOCK AND STOCKING AID

The Brevetti Sock and stocking aid is designed to help those with restricted mobility to apply socks and stockings easily. The Brevetti consists of a coated steel tube with a plastic collar to hold the garment open on one end and a plastic shoehorn with a notch for easy sock removal on the other.

Specifications:

Length	73.7 cm (29")
Weight	386 g (13.6 oz)



Code	Description
091095603	Brevetti Sock and Stocking Aid



DORKING STOCKING DONNER

The Dorking Stocking Donner is a compact, durable tool that extends the user's reach to help them dress independently. The extendable, hold-down handles rotate 90°, creating a comfortable position for the user while applying stockings. The central column, which holds the stocking, has a small ridge that prevents movement and widens at the top to accommodate the larger calf. The epoxy coating creates a smooth surface, reducing the risk of snagging.

Specifications:

Height	31 cm (12¼")
Width	20.5 cm (8")
Column height	20.5 cm (8")
Column width	12.5 cm (5")
Column depth	8.5 cm (3½")
Height with handles	56 cm (22")
Weight	1.1 kg (2½ lb)



Code	Description
081330869	Dorking Stocking Donner



DORKING TIGHTS DONNER

The Dorking Tights Donner is like the Dorking Stocking Donner but with twin columns designed for use with tights. The lightweight frame enables those who need support with upper body strength to apply the tights easily.

Specifications:

Width	38 cm (15")
Height with handle	56 cm (22")

Code	Description
091095447	Dorking Tights Donner



CONTOUR TURNER

The Contour Turner offers leverage for turning small household items, such as gas or radiator taps. The head contains spring-loaded stainless-steel rods, which when pressed around an object, will retract and conform around it, creating a grip around the object when the handle is turned. The handle is made of easy-to-grip moulded plastic.

Specifications:

Head size — width	9.5 cm (3¾")
Head size — depth	5 cm (2")
Weight	130 g (4.6 oz)

Code	Description
081499623	Counter Turner



STEP STOOL

The Step Stool is a simple, sturdy stool consisting of a chrome-plated steel frame with non-slip rubber mats. Rubber ferrules provide stability on floor surfaces. Available with or without a rail.

Specifications:

Top area	28 x 36 cm (11 x 14¼")
Base area	30.5 x 41 cm (12 x 16¼")
Base height	22.5 cm (9")
Maximum load	120 kg (18¾ st; 264½ lb)

Code	Description
091102680	Step Stool
091102698	Step Stool with Handrail



PLUG PULL

The Plug Pull is a self-adhesive, plastic moulded T-piece that fixes to the back of a standard, three-pin plug, providing leverage to remove the plug from the socket. It is also helpful for drawers and doors where extra leverage is required. Supplied in a pack of three.

Code	Description
091102441	Plug Pull — Set of 3



PEN AND PENCIL HOLDER

The Pen and Pencil Holder is constructed from soft PVC moulding, which slides onto pens and pencils to make gripping and controlling writing easier. Supplied in a pack of three.

Specifications:

Length	73.7 cm (29")
Weight	386 g (13.6 oz)

Code	Description
091104603	Pen and Pencil Holder — Set/3



MEDI-DISPENSER

The Medi-Dispenser is an organizer that holds a week's medication with as many as four daily doses. The plastic wallet has seven plastic containers printed with the day of the week and embossed with the number of the day of the week.

Code	Description
091146398	Medi-Dispenser



REVOREACH™ EXCEL

The RevoReach™ Excel is a robust reacher with an ergonomically designed handle and trigger that require minimal finger extension, allowing the use of all fingers for a powerful grip. The short trigger is ideal for those with arthritic fingers or a weak grip.



Revolving jaw head prevents the need to rotate wrist



Slip resistant fingerprint-lined jaws for optimum grip



Code	Description	Length	Weight
091204825	Standard	66 cm (26")	210 g (7.4 oz)
091204833	Long	81 cm (32")	230 g (8.1 oz)



REVOREACH™ GRIPLOCK

The RevoReach™ Griplock is like the RevoReach™ Excel and helps those who need support with dexterity. The Griplock design features a locking mechanism that allows the user to grip and hold an object and then engage the lock after they have achieved a secure hold, automatically locking as the user squeezes the trigger.



Code	Description	Length	Weight
091204791	Standard	66 cm (26")	220 g (7.8 oz)
091204809	Long	81 cm (32")	240 g (8.5 oz)



HANDI-REACHER™

The Handi-Reacher™ is an ultra-lightweight and robust reacher that assists those with impaired reach and mobility. Featuring a single moving jaw operated with a trigger action from the handle, the Handi-Reacher™ comes with a magnet at the end for picking up small metal objects. The jaw rotates 360° to avoid the need to turn the wrist. The trigger — operated by four fingers instead of the typical two — is easy to grip and use.



Swivelling,
rubber-lined
jaw



Full hand-
operated
trigger



Code	Description	Length	Weight
091114420	Extra Long	90 cm (35½")	220 g (7.8 oz)



PICKUP REACHER

The Pickup Reach is a lightweight, robust reacher ideal for retrieving or reaching various items in the home, garden or stores. Available in four lengths and with optional extras, it has twin-gripping jaws with sculptured liners and a trigger designed for comfort. The pickup is operated by a wire, not a cord, eliminating fraying and stretching. A magnetic head enables easy retrieval of pins, paper clips and other small, metal items. The handle is large and roomy to allow for the use of all fingers, and the short trigger distance is ideal for those with stiff, arthritic fingers or limited mobility. It has an easily rotatable head and a removable stick clip to connect the reacher to wheelchairs, walking frames and other mobility devices.

Specifications:

	Length	Weight
Standard	60 cm (24")	215 g (7.6 oz)
Long	75 cm (30")	235 g (8.3 oz)
Extra long	90 cm (35½")	250 g (8.8 oz)



Code	Description
081306430	Standard
091114511	Long with Forearm Support
081306497	Extra Long
091114503	Long Retail Packed



Bed & Seating equipment



OVER ARMCHAIR TABLE

The Over Armchair Table has a curved, laminated teak tabletop that is easy to adjust in height and angle. The modern style includes a small table that can be fitted to the main table's left or right side.

Specifications:

Table top	109 x 41 cm (43 x 16")
Side table top	20 x 20 cm (8 x 8")
Overall width	13.1 cm (51½")
Height adjustment	75-91.5 cm (29½-36")
Internal width	90.6 cm (35½")
Maximum load	10 kg (22 lb)



Code	Description
081610971	Over Armchair Table



LAMINATED TOP OVERBED TABLE

The Laminated Top Overbed Table has a teak top and two independent points to adjust the height to angle the top. Tabletop rails prevent items from slipping off. The table does not have castors, so it can be used with a bed or chair with low clearance. Available fixed or flat-packed.

Specifications:

Height adjustment	61-91 cm (24-36")
Table top size	61 x 40 cm (24 x 15¾")
Maximum load	10 kg (22 lb)

Code	Description
081611037	Fixed
091558253	Flat-packed



081611037

091558253



LAMINATED TOP OVERBED TABLE WITH CASTORS

The Laminated Top Overbed Table with Castors is like the Laminated Top Overbed Table but with castors to easily move the table as needed. Available fixed or flat-packed.

Specifications:

Height adjustment	61-91 cm (24-36")
Table size	61 x 40 cm (24 x 15¾")
Castors	5 cm (2")
Maximum load	10 kg (22 lb)

Code	Description
091325745	Fixed
091558261	Flat-packed



PLASTIC TOP OVERBED TABLE

The Plastic Top Overbed Table features a lightweight plastic top for easy cleaning and care. A ridge around the table edges helps prevent items from slipping off and can catch small spills. The ratchet handle adjusts the height, and the table can be easily moved with four castors.

Specifications:

Height adjustment	75-106 cm (29½-41¾")
Tabletop	70.5 x 37.5 cm (27¾ x 14¾")
Clearance required	9.25 cm (3¾")
Width between bars — base	56 cm (22")
Width — base bar	34 cm (13⅓")

Code	Description
081611045	Plastic Top Overbed Table



LAMINATED TEAK TOP OVERBED TABLE

The Laminated Teak Top Overbed Table is manufactured with a chrome-plated and powder-coated steel tubing frame and features an easy-to-clean teak top. Two independent height adjustments allow the top to be angled, while the raised edges prevent items from slipping off and can catch small spills. The table can be easily moved with four castors.

Specifications:

Height adjustment	53-82 cm (21¼-32¾")
Tabletop size	54 x 40 cm (21¼ x 15¾")
Maximum load	10 kg (22 lb)

Code	Description
091439157	Laminated Teak Top Overbed Table



ADJUSTABLE WOODEN BED TRAY

The Adjustable Bed Tray is an attractive, practical wooden-frame bed tray that can be adjusted to five angles for reading, writing and other activities while sitting in bed. The legs, which are laminated for easy cleaning, fold out to allow the tray to fit comfortably across the lap.

Specifications:

Height	23.5 cm (9¼")
Width	54 cm (21¼")
Depth	31 cm (12¼")



Code	Description
091100361	Adjustable Wooden Bed Tray

WALKER TROLLEY

The Walker Trolley is a sturdy walking aid that helps the user transport items around the home. The easy-to-clean design includes clip-on plastic shelves and lockable hand brakes. Rails on the tray edges help prevent objects from slipping off.

Specifications:

Overall width	58.5 cm (23")
Overall depth	56 cm (22")
Handle height	75.5-90.5 cm (29¾-35¾")
Wheel diameter	15 cm (6")
Maximum load	10 kg (22 lb)



091557289



081449230

Code	Description	Colour
091557289	Walker Trolley	White
081449230	Walker Trolley	Beige



HEIGHT-ADJUSTED TROLLEY

Height-adjustable Trolley features a metal frame with plastic clip-on shelves measuring 46 x 33 cm (18 x 13") that are easy to remove for cleaning. The trolley is height-adjustable and has four castors. Supplied flat-packed or fixed.

Specifications:

Height adjustment	88-100 cm (34½-39½")
Maximum load	8 kg (17¾ lb)

Code	Description
091438969	Flat-packed
091435213	Fully assembled



WOODEN CHAIR RAISERS

Wooden Chair Raisers are cubes in bamboo material with three 1.6 cm (5/8") recesses of different diameters to suit various leg widths. Supplied in a set of four.

Specifications:

Height	9.3 cm (3¾")
Height raise	7.5 cm (3")
Maximum diameter of leg	3, 5 or 7.5 cm (1¼, 2 or 3")
Weight	425 g (15 oz) each

Code	Description
091088855	Wooden Chair Raisers — Set of 4



STACKABLE CONE RAISERS — BED OR CHAIR

Stackable Cone Raisers are lightweight but strong moulded cones with a broad base for stability. Each features a recess for the chair or bed leg and anti-slip pads to hold the raiser in place. Set of four available for bed or chair.

Specifications:

	Chair	Bed
Height	9.3 cm (3¾")	14 cm (5½")
Height raise	7.5 cm (3")	7.5 cm (3")
Maximum diameter of leg	3, 5 or 7.5 cm (1¼, 2 or 3")	3, 5 or 7.5 cm (1¼, 2 or 3")
Weight	105 g (3.7 oz)	145 g (5.2 oz)
Maximum load	222 kg (35 st; 489½ lb)	222 kg (35 st; 489½ lb)

Code	Description
091088814	Stackable Cone Chair Raisers — Set of 4
091088830	Stackable Cone Bed Raisers — Set of 4



LOOP BED GRAB RAIL

The Loop Bed Grab Rail provides extra support when getting into and out of bed and is an aid in achieving a comfortable position. The rail is positioned onto the bed frame and secured with two ergonomic, easy-to-grip hand knobs.

Specifications:

Gap at widest point 9.8 cm (3⁷/₈")

Code	Description
091439546	Loop Bed Grab Rail



DOUBLE LOOP BED GRAB RAIL

The Double Loop Bed Grab Rail is designed with a loop on each side of the mattress to assist users when getting into and out of bed. It fits securely under the mattress and is width-adjustable to fit single and double beds.

Specifications:

Length 98-143 cm (38¹/₂-56¹/₄")

Height 44.5 cm (17¹/₂")

Code	Description
091325661	Double Loop Bed Grab Rail



SLATTED BED GRAB RAIL

The Slatted Bed Grab Rail is suitable for slatted beds, with a large metal base frame for maximum under-mattress stability.

Specifications:

Base width 62.5 cm (24¹/₂")

Base length 77.5 cm (30¹/₂")

Base to top of handle 48 cm (19")

Code	Description
091325752	Slatted Bed Grab Rail



BED GRAB RAIL WITH LEGS

The Bed Grab Rail with Legs is ultra-sturdy, offering stability and confidence to those who have difficulty getting into and out of bed. The height-adjustable frame can be positioned to fit any divan or similar bed style. Two legs rest on the floor, and three colour-coded straps (supplied) secure the frame. Easy to assemble.

Specifications:

Handle height	85-95 cm (33½-37½") from the floor
Handle width at top	26.5 cm (10½")
Size of under-mattress frame	62 x 32.5 cm (24½ x 12¾")
Weight	3.5 kg (7¾ lb)



Code	Description
091089481	Bed Grab Rail with Legs



COT SIDE BUMPERS

Cot Side Bumpers are made from thick foam and covered with wipe-clean fabric. They are made to fit various cot sizes and easily fasten into place. Supplied in a set of two.

Specifications:

Length	152.4 cm (60")
Depth	38.1 cm (15")
Weight	1 kg (2¼ lb) each



Code	Description
091089333	Bed Rail Buffer Pads — Set of 2



ADJUSTABLE BACK REST

The Adjustable Back Rest is secure and comfortable and will adjust to many positions, from low-lying angles to upright sitting positions. Breathable nylon fabric covers the metal frame. Includes a small pillow for additional comfort.

Specifications:

Maximum height of back	58 cm (22¾")
Width	65 cm (25½")
Depth	50 cm (19½")
Weight	2.8 kg (6¼ lb)
Maximum load	113 kg (17¾ st; 249¾ lb)

Code	Description
091089721	Adjustable Back Rest



CHROME FOLDING BED CRADLE

The Chrome Folding Bed Cradle supports the weight of bedclothes to reduce contact with the user. It is manufactured from chrome-plated steel and designed to fit under the mattress and fold away flat for storage.

Specifications:

Height	51 cm (20")
Length	50 cm (20½")
Overall width	21 cm (8¼")
Maximum load	5 kg (11 lb)

Code	Description
091325620	Chrome Folding Bed Cradle



RIGID LEG LIFTER

The Rigid Leg Lifter is a webbing strap with a stiffened stem that enables the user to move an immobile leg when in a bed or wheelchair or using a footstool. The upper loop fits around the hand and wrist. The lower loop slides over the foot to offer control in lifting the leg.

Specifications:

Length	86.4 cm (34")
Weight	78 g (2.7 oz)

Code	Description
091118793	Rigid Leg Lifter



STIRRUP-STYLE LEG LIFTER

The Stirrup-style Leg Lifter is similar in style and function to the Rigid Leg Lifter but has a stiffened triangular stirrup, making it easier to fit over shoes and offering more foot support.

Specifications:

Length 91.5 cm (36")

Weight 120 g (4.2 oz)

Code	Description
091118801	Stirrup-style Leg Lifter



DUAL-HANDLE LEG LIFTER

The Dual Handle Leg Lifter is designed to help users reposition their legs when they have limited or no use of their legs. Two handles enable users to utilise both arms when moving the leg, and the sturdy centre rod stabilises the device. Ideal for using while in a chair or bed.

Specifications:

Length 86.5 cm (34")

Weight 110 g (3.9 oz)

Code	Description
091118819	Dual-Handle Leg Lifter



ROPE LADDER HOIST

The Rope Ladder Hoist offers a simple way to assist with sitting up in bed. The bottom ends of the cords tie around the bed legs, and plastic ladder rungs enable the user to pull themselves up to a sitting position with a hand-over-hand movement. Suitable for single or double beds.

Specifications:

Rung length 20 cm (8")

Diameter 2.5 cm (1")

Cord length 3 m (118")

Weight 285 g (10 oz)

Maximum load 60 kg (9½ st; 132¾ lb)

Code	Description
091089895	Rope Ladder Hoist



PLASTIC CURVED TRANSFER BOARD

The Plastic Curved Transfer Board enables various positions to assist with functional transfers for individuals with limited strength and mobility. The transfer board is designed to aid seated transfers of individuals between two surfaces, such as bed to wheelchair or wheelchair to chair or toilet. The board helps reduce the risk of injury and strain on the caregiver assisting with transfers.

Specifications:

Length	71 cm (27½")
Width	37 cm (14½")
Weight	1.95 kg (4¾ lb)



Code	Description
081306554	Plastic Curved Transfer Board

152 kg
23¾ st 335¾ lb
Maximum user weight

WOODEN TRANSFER BOARD

The durable Wooden Transfer Board is ideal for bed, wheelchair and vehicle transfers. The top layer is varnished to facilitate a friction-free transfer, reducing strain and eliminating the need to lift patients. The chamfered ends offer secure positioning during transfers.



091120633



091120724



Code	Length
091120633	60 cm (24") without hand holes
091120724	76 cm (30") with hand holes

Weight
1.1 kg (2½ lb)
1.38 kg (3 lb)

190 kg
30 st 419 lb
Maximum user weight



ORBI-TURN TRANSFER AID

The Orbi-Turn Transfer Aid features a steel frame and is designed for caregivers to facilitate standing and seat-to-seat transfers while minimising the risk of injury to the patient and themselves. The Orbi-Turn enables a patient with moderate upper and lower limb strength to participate in functional transfers while ensuring good working postures for caregivers. A height-adjustable, angled knee pad facilitates typical biomechanical hip and knee movement.

Specifications:

	Lowest Position	Highest Position
Bottom handle	81.5 cm (32")	117.5 cm (46¼")
Top handle	92 cm (36¾")	131.5 cm (51¾")
Knee pad (mid point)	23.5 cm (9¼ lb)	51 cm (20")



Code	Description
081590769	Orbi-Turn Transfer Aid

190 kg

30 st 419 lb

Maximum user weight



Kitchen & Dining



QUEENS CUTLERY



Queens Cutlery includes a range of stainless-steel utensils with lightweight, ABS plastic built-up handles designed to help those with weak or restricted grips. Each handle is cylindrical with a slight taper away from the utensil and smooth contours.

- **Queens Standard Cutlery** is a popular range of standard utensils mounted in the Queens built-up handle. The knife, fork, spoon, junior spoon and teaspoon are available individually or as a set. Dishwasher-safe.
- **Queens Angled Cutlery** features stainless steel utensils with angled heads to help individuals with limited wrist or arm movement. The fork and spoon are available as left- or right-hand options. Dishwasher-safe.
- **One-handed Cutlery** features a stainless steel utensil with the Queens built-up handle and combines multiple cutlery functions. They are designed to help amputees and people with semi-paralysis after a stroke eat independently.

Specifications:

Handle length	9.8 cm (4 $\frac{5}{8}$ ")
Handle diameter	3 cm (1 $\frac{1}{4}$ ")

Code	Description	Weight
Standard Cutlery		
081305820	A — Knife	45 g (1.6 oz)
081305838	B — Fork	45 g (1.5 oz)
081305846	C — Spoon	45 g (1.5 oz)
081305853	D — Junior Spoon	45 g (1.5 oz)
091097427	E — Teaspoon	45 g (1.5 oz)
091097799	Queens Retail Packed Set including Knife, Fork, Spoon and Junior Spoon (does not include teaspoon)	
Angled Cutlery		
081305879	Knife (not illustrated)	45 g (1.6 oz)
081305895	F — Right-handed Fork	55 g (1.9 oz)
081305887	G — Right-handed Fork Ivory	55 g (1.9 oz)
081305911	H — Right-handed Spoon	55 g (1.9 oz)
081305903	I — Left-handed Spoon	55 g (1.9 oz)
One-handed Cutlery		
081305861	J — Spork	45 g (1.5 oz)
091097468	K — Nelson Knife	55 g (1.9 oz)




CARING CUTLERY

Caring Cutlery is a set of stainless steel utensils featuring contoured plastic handles with built-up ends. These utensils are ideal for individuals with weak grips and restricted wrist or finger movement. The knives and forks have a shaped indent on the top for the index finger to help with directional control. Also available as a four-piece set including knife, fork, spoon and teaspoon.

Specifications:

Handle length	12.7 cm (5")
Utensil weight	45-60 g (1.6-2.1 oz)

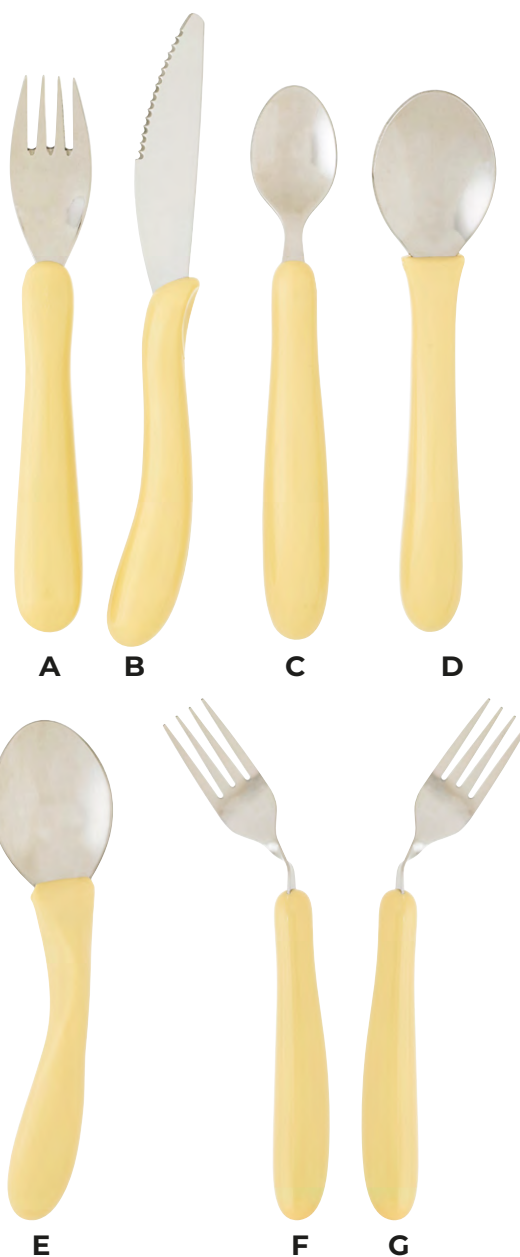

Code	Description	Weight
Standard Cutlery		
091098607	A — Knife, Ivory	50 g (1.8 oz)
091098623	B — Fork, Ivory	50 g (1.8 oz)
091098631	C — Spoon, Ivory	60 g (2.1 oz)
091098649	D — Teaspoon, Ivory	50 g (1.8 oz)

Angled Cutlery

091098672	E — Right-handed Spoon, Ivory	60 g (2.1 oz)
091098680	F — Left-handed Fork, Ivory	50 g (1.8 oz)
091098698	G — Right-handed Fork Ivory	50 g (1.8 oz)

Cutlery Sets

081533801	Ivory Set including Knife, Fork, Spoon and Teaspoon	215 g (7.6 oz)
081305994	Black Set including Knife, Fork, Spoon and Teaspoon	215 g (7.6 oz)




LIGHTWEIGHT FOAM-HANDLED CUTLERY

Lightweight Foam-handled Cutlery offers a comfortable design to help those with weak or restricted grips. The stainless steel utensils feature light grey foam handles that are ribbed and soft. Hand-washable only.

- **Standard Cutlery** including a knife, fork, spoon and teaspoon, is available individually or as a set. Hand-washable only.
- **Angled Cutlery** features stainless steel utensils with angled heads to help individuals with limited wrist or arm movement. The fork and spoon are at 40° angles and available as left- or right-hand options. The knife is at a 90° angle and can be used without cocking the wrist.
- **One-handed Cutlery** features a stainless steel utensil with the lightweight foam handle that combines multiple cutlery functions. They are designed to help amputees and people with semi-paralysis after a stroke eat independently.

Specifications:

Handle length	11.5 cm (4½")
Handle diameter	2.8 cm (1¼")

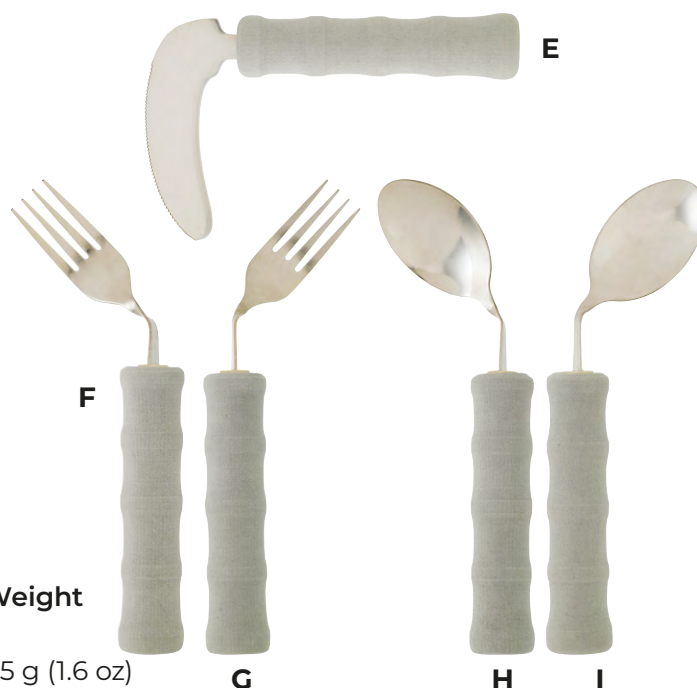

Code	Description	Weight
Standard Cutlery		
081305986	A — Knife	45 g (1.6 oz)
081305994	B — Fork	34 g (1.2 oz)
081306000	C — Spoon	42 g (1.5 oz)
081306018	D — Teaspoon	30 g (1.1 oz)
091098516	Packed Set — 1 of each	151 g (5.3 oz)

Angled Cutlery

081306034	E — Right-handed Spoon, Ivory	45 g (1.6 oz)
091098540	F — Left-handed Fork, Ivory	40 g (1.4 oz)
091098532	G — Right-handed Fork Ivory	40 g (1.4 oz)
091098565	H — Right-handed Spoon	46 g (1.6 oz)
091098557	I — Left-handed Spoon	46 g (1.6 oz)

One-handed Cutlery

091098466	J — Splayed Fork	40 g (1.4 oz)
091098441	K — Nelson Knife	45 g (1.6 oz)



NEWSTEAD CUTLERY

Newstead Cutlery features ergonomic, flexible, moulded plastic handles designed to adapt to the user's grip. With a comfortable rounded shape, it fits naturally into the palm and is ideal for people with arthritis, neurological impairments and other conditions that cause poor grip. All pieces, made of stainless steel and containing no latex, are dishwasher safe.

- **Standard cutlery set includes a knife, fork, spoon and teaspoon.**
- **Angled cutlery features angled heads to help individuals with limited wrist or arm movement. The fork and spoon are at 40° angles, available as left- or right-hand options and may be bent to adjust the angle. The angled knife can be used for cutting without cocking the wrist.**
- **Weighted cutlery features standard cutlery heads but with a heavier handle, which helps reduce tremors to reduce spills.**
- **One-handed cutlery make dining with one hand easier for people with limited hand movement after a stroke, amputees and others who need to eat with one hand.**

Specifications:

Standard

Weight	80 g (2.8 oz) per utensil
Set	320 g (11.3 oz)

Weighted

Weight	130 g (4.6 oz) per utensil
Set	520 g (1¼ lb)

Code	Description	Weight
Standard Cutlery		
091097872	A — Cutlery Set	320 g (11.3 oz)

Angled Cutlery

091097880	B — Right-handed Fork	80 g (2.8 oz)
091097898	C — Right-handed Spoon	80 g (2.8 oz)
091097906	D — Left-handed Fork	80 g (2.8 oz)
091097914	E — Left-handed Spoon	80 g (2.8 oz)
091097922	F — Angled Knife	80 g (2.8 oz)

Weighted Cutlery


091097948	G — Weighted Knife	130 g (4.6 oz)
091097955	H — Weighted Fork	130 g (4.6 oz)
091097963	I — Weighted Spoon	130 g (4.6 oz)
091097971	J — Weighted Teaspoon	130 g (4.6 oz)
091097989	Weighted set	520 g (18.3 oz)

One-handed Cutlery

091097997	K — Splayed Fork	80 g (2.8 oz)
091098003	L — Nelson Knife	80 g (2.8 oz)



A — complete set including Knife, Fork, Spoon and Teaspoon

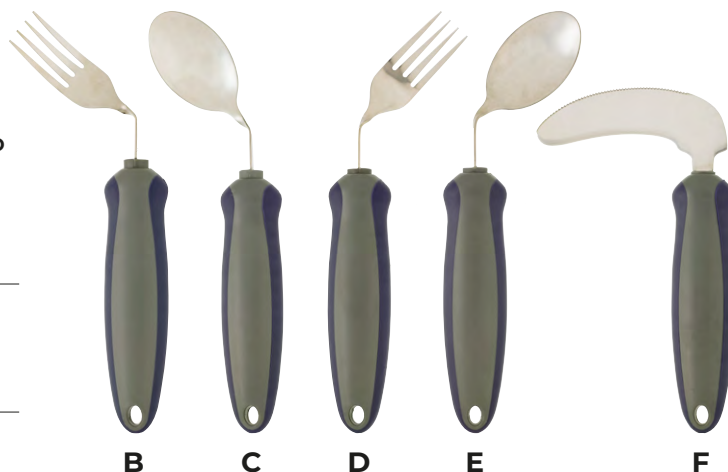

CARING MUG

The Caring Mug has a wide, stable base with two large contoured and angled handles, which enable it to be held in the correct position for drinking. The lid has a uniquely designed spout and matches the Caring Cutlery range. It is useful for those with tremors or who cannot sit upright. The mug is microwave- and dishwasher-safe. Supplied with a small aperture for liquids or a large aperture for semi solids.

Specifications:

Capacity 300 ml (10¼ fl oz)

Code	Description
091099555	Small aperture
091165422	Large aperture



HANDLED MUG WITH LIDS

The Handled Mug with Lids is a popular mug with large handles and shaped finger holds for ease of use. Liquid consumption is easily measured by marked gradations. The mug is dishwasher-safe, but the lids should be cleaned in warm, soapy water. Supplied with one spout lid and one splash lid.

Specifications:

Capacity 270 ml (9 fl oz)

Code	Description
091099571	Handled Mug with Lids — Set of 2



FEEDING CUP

The Feeding Cup is designed for easy gripping and features a vented lid with a drinking spout. The cups are dishwasher-safe up to 135° C (275° F). Supplied in pairs.

Specifications:

Capacity 250 ml (9 fl oz)
Weight 80 g (3 oz)

Code	Description
091099977	4 mm aperture
091099985	8 mm aperture



SHATTERPROOF MUG

The 400-ml (13½ fl oz) Shatterproof Mug is made from a tough material to prevent accidental breakage. Liquid consumption is easily measured by marked gradations, and the mug is supplied with a choice of lids. The spout lid helps with drinking, and the recessed, spill-proof lid avoids splashes. The mug is microwave- and dishwasher-safe, but the lids should be hand-washed.



Code	Description
081569656	Shatterproof Mug

WIDE-BASE MUG WITH WHITE SPOUT LID

The Wide-base Mug with White Spout Lid is made from high-strength polycarbonate and has a wide base for stability to help prevent tipping. The dual handles provide a secure grip, making it ideal for anyone who has a limited grasp. Liquid consumption is easily measured by marked gradations. It is dishwasher-safe. Also available with a set of three transparent lids.

Specifications:

Capacity	285 ml (9½ fl oz)
----------	-------------------



Code	Description
091100049	Mug with Solid White Lid
091001262	Mug with Transparent Lid



INCURVE PLATE SURROUND

The Incurve Plate Surround is a flexible guard that can be fitted to plates to assist with one-handed eating. Three clips attach it to the rim and create a high inward sloping face to help with food collection and prevent spills. It is microwave- and dishwasher-safe.

Specifications:

External diameter	19-24.5 cm (7½-10")
Weight	45 g (1.6 oz)



Code	Description
081306091	Incurve Plate Surround




REFLEX KITCHEN KNIFE COLLECTION

Reflex Kitchen knives offer a secure but easy-to-grip handle, which is contoured and comfortable. The handle has been designed to counterbalance the weight of the knife. The angled, stainless steel blade minimises the user's wrist movement, making this collection ideal for those with weakened grip strength. Dishwasher-safe.

Specifications:

	Blade length	Weight
Reflex Carving Knife	20.3 cm (8")	155 g (5.5 oz)
Reflex Chef's Knife	20.3 cm (8")	165 g (5.8 oz)
Reflex Preparation Knife	12.7 cm (5")	140 g (4.9 oz)
Reflex Slicing Knife	25.4 cm (10")	125 g (4.4 oz)

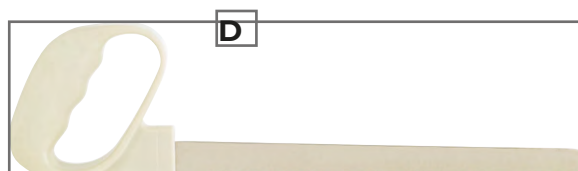

Code	Description
081305788	A — Carving Knife
081305796	B — Chef's Knife
081305770	C — Preparation Knife
081305804	D — Slicing Knife



LONG-HANDLED BOTTLE HOLDER

The Long-handled Bottle Holder is a wire frame that helps users carry bottles quickly and safely without bending to pick them up. This leaves one hand free for support.

Specifications:

Handle length	60 cm (24")
Max bottle diameter	8.3 cm (3¼")
Weight	610 g (1½ lb)

Code	Description
091096916	Long-handled Bottle Holder



UNIVERSAL KETTLE TIPPER

The Universal Kettle Tipper features a classic design compatible with most kettle types, including standard, jug and cordless. The lever can be positioned on either side to tip the kettle without having to lift it. For added safety, the kettle rocks back into an upright position when released. The base plate has a front sliding stop so the user can position the kettle where it is needed. It is then secured with two sturdy straps. The tipper is manufactured from attractive, heat-resistant plastic for durability and has rubber feet to prevent sliding on the work surface.

Specifications:

Base plate diameter	20 cm (7⅞")
Weight	190 g (6.7 oz)



Code	Description
091207281	Universal Kettle Tipper



PLASTIC SPREADING AND CHOPPING BOARD

The Spreading and Chopping Board is a durable, one-piece, moulded board with one raised, L-shaped corner, allowing jams and other foods to be spread on bread using one hand. It features six stainless steel spikes that can be inserted into the board and are designed to hold food while slicing or chopping.

Specifications:

Size 18.1 x 25.6 cm (7 $\frac{1}{8}$ x 10 $\frac{3}{8}$ ")

Weight 190 g (6.7 oz)



Code

081569623

Description

Plastic Spreading and Chopping Board


PLASTIC SPREADING BOARD

The Plastic Spreading Board is a durable, one-piece moulded board with one raised, L-shaped corner, assisting the user when spreading jams and other foods on bread using one hand.

Specifications:

Size 18.1 x 25.6 cm (7 $\frac{1}{8}$ x 10 $\frac{3}{8}$ ")

Weight 190 g (6.7 oz)

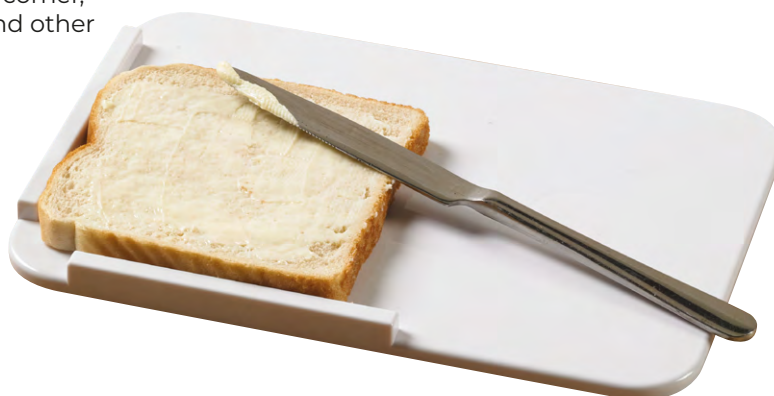

Code

091096742

Description

Plastic Spreading Board

STAINLESS STEEL COOKING BASKET

The Stainless Steel Cooking Basket enables users to cook vegetables without ever needing to lift a pan full of boiling water. Place vegetables into the basket, insert them into the pan, cook the food and remove the basket when complete. The large, flame-retardant, moulded nylon handle is easy to grip and offers heat insulation.

Specifications:

Top diameter 16.5 cm (6 $\frac{1}{2}$ ")

Bottom diameter 12 cm (4 $\frac{1}{2}$ ")

Depth 9 cm (3 $\frac{1}{2}$ ")

Weight 245 g (8.6 oz)



Code

091096890

Description

Stainless Steel Cooking Basket



KITCHEN WORKSTATION



The Kitchen Workstation is a multifunction food preparation unit that assists with many day-to-day kitchen tasks. The large, durable chopping board features a raised, L-shaped corner that can hold bread and toast in place while spreading, and removable stainless steel spikes secure bread, fruit or vegetables while being sliced, chopped or peeled. A removable clamp can hold various objects in place, such as a bowl, a tin or food, making food preparation easier. The removable grater and slicer mount to the

chopping board, making one-handed use possible with minimum mess and waste. Four non-slip suction pads firmly secure the workstation in place. Dishwasher-safe.

Specifications:

Overall dimension	50 x 30 cm (19¾ x 11¾")
Clamp width	11 cm (4¼")
Weight	1.24 kg (2¾ lb)

Code	Description
081112259	Kitchen Workstation

EVERYDAY BIB

The Everyday Bib is practical for everyday use. The colourfast polycotton outside is durable but absorbent. The soft and flexible waterproof backing helps reduce clothes stains. All bibs are fitted with press studs, so the user can form a pocket to catch excess food. Users can fasten the bibs with hook-and-loop patches for quick, easy donning or metal poppers for extra security, with two alternative fastening positions for the poppers. The bibs are fully washable and can be tumble-dried at low temperatures. Available in three sizes and two colours.

Specifications:

Size	
Small	58 cm x 38 cm (23" x 15")
Medium	68 cm x 46 cm (26¾" x 18")
Large	89 cm x 46 cm (35" x 18")

Code	Description
091421148	Everyday Bib, Red, Small
091421155	Everyday Bib, Blue, Small
091099399	Everyday Bib, Blue, Medium
091099407	Everyday Bib, Red, Medium
091099365	Everyday Bib, Blue, Large
091099373	Everyday Bib, Red, Large



FREEHAND TRAY

The Freehand Tray is designed to enable one-handed carrying, leaving a hand to open doors, climb stairs and provide support. The unique handle folds flat inside the inner rim of the tray and lifts and locks into place in the upright position for carrying. The handle is easily removed for cleaning or for use as a standard tray. Dishwasher safe up to 130°C (266°F).

Specifications:

Size	46.5 cm x 30.5 cm
Weight	1.1 kg (2½ lb)
Max load	4 kg (8¾ lb)



Code	Description
091100320	Freehand Tray with Carry Handle

STAY TRAY WITH BEAN BAG

The Stay Tray with Bean Bag sits comfortably on your lap to provide a secure surface for eating, drinking or recreational activities, making it ideal for use on a chair or sitting in bed. The durable plastic tray has built-in handles and comes with a detachable bean bag base, which is removable and can be machine-washed. Includes a non-slip, heat-resistant silicon mat to stop plates or cups from sliding.

Specifications:

Size	18.1 x 25.6 cm (7¼ x 10¼")
Weight	190 g (6.7 oz)



Code	Description
091100338	Stay Tray with Bean Bag



Independence in the Home

Homecraft® was founded in the 1940s to honour soldiers and armed forces veterans, offering dignity to those who need extra support on life's journey. Although times have changed, the fundamental values of Homecraft have remained steadfast. Today, Homecraft fosters independence for individuals of all ages with an extensive selection of daily living aids available globally, including bathing, toileting, kitchen aids and more. Whether you are a clinician, retailer, consumer or family, Homecraft provides peace of mind to live life to the fullest.

Be your ideal self

Homecraft is at the heart of every home when you need us. Our mission is to empower individuals in their environments. Whether it is buttoning a shirt, taking a bath or using the toilet, our aids to daily living (ADLs) are there to support you.

We know daily tasks can be challenging. Homecraft was built with the spirit and belief that each carefully crafted product will empower you to lead your own life. People turn a house into a home, and Homecraft embraces your safe space with designs to fit your lifestyle.





Performance Health

INTERNATIONAL SALES OFFICE

208 Kifissias Avenue, 15231 Halandri - Athens - Greece

tel.: +30 210 6849 970, +30 210 6712 266

fax: +30 210 6829 255

info@performancehealth.com

CORPORATE HEADQUARTERS

1245 Home Avenue, Akron - Ohio 44310 - USA

tel.: +1 330 633 8460

